

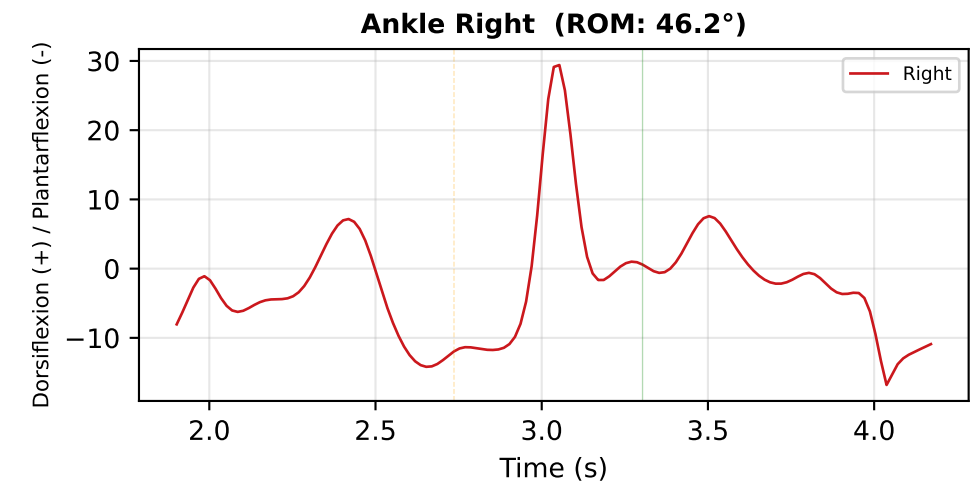
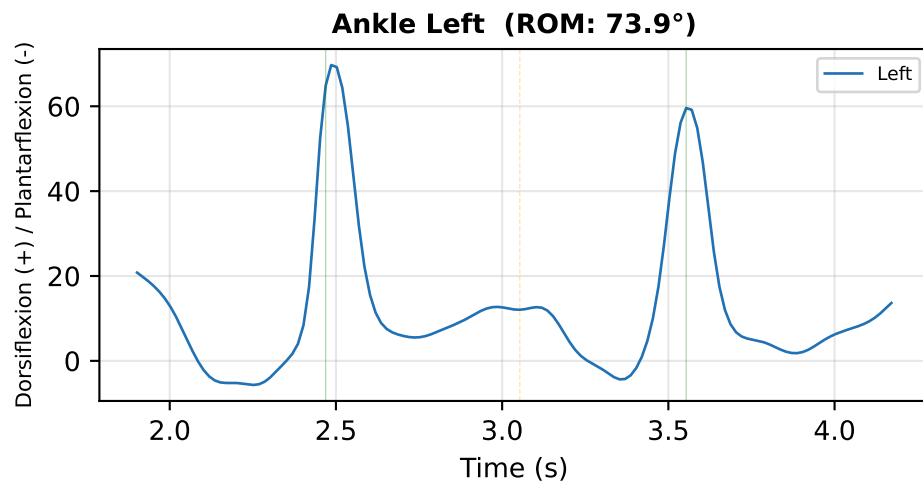
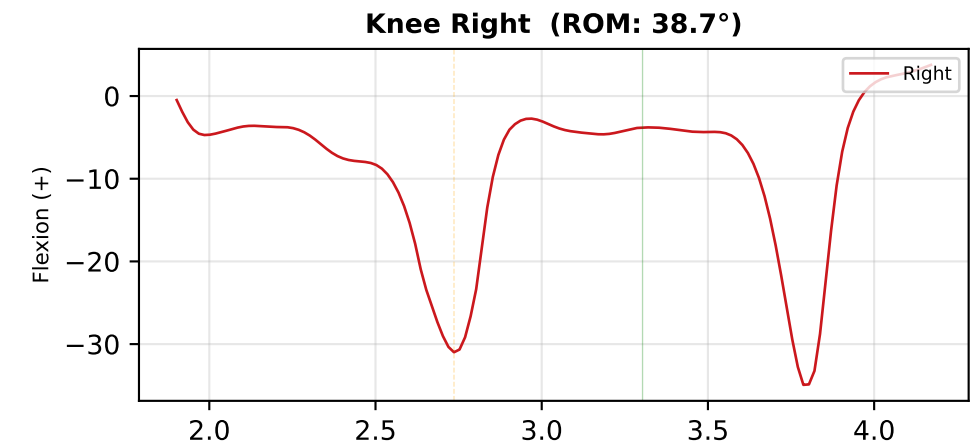
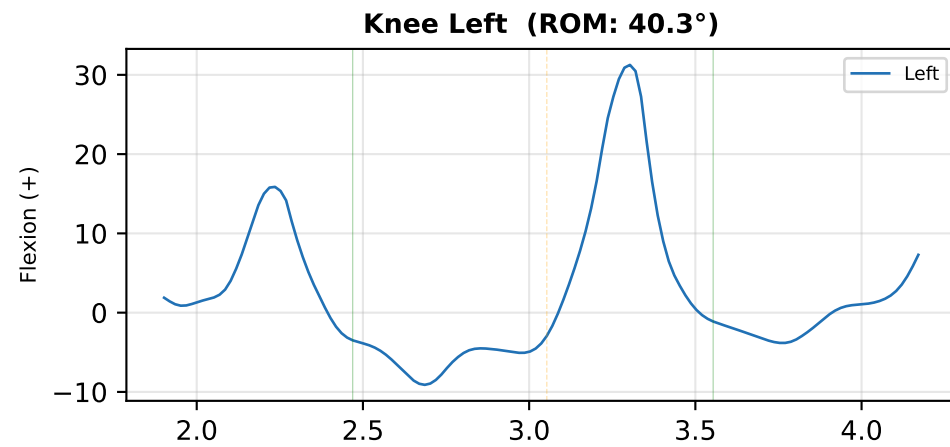
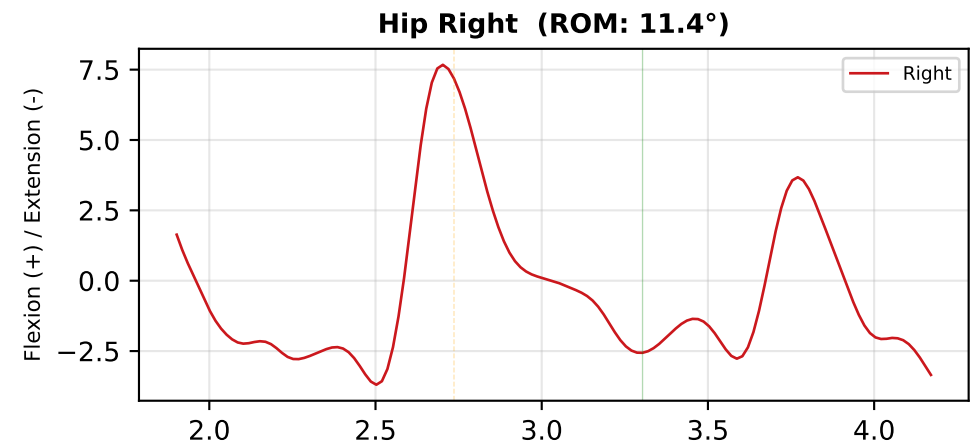
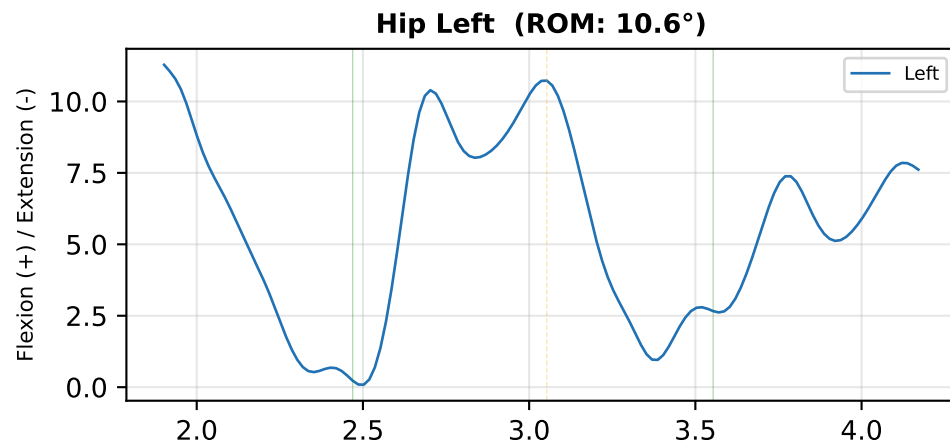
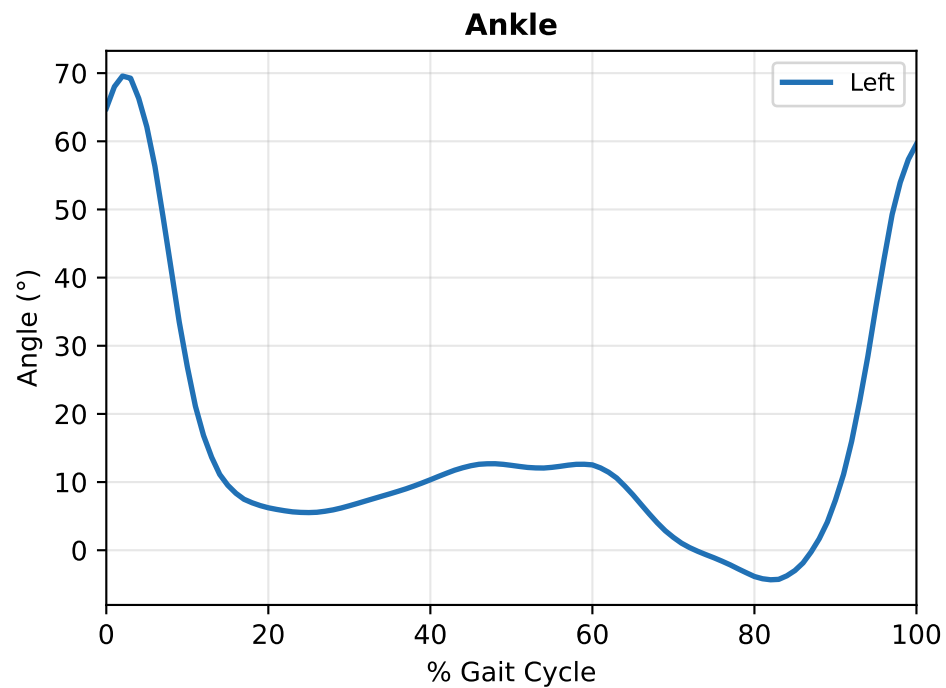


Overview — Joint Angles

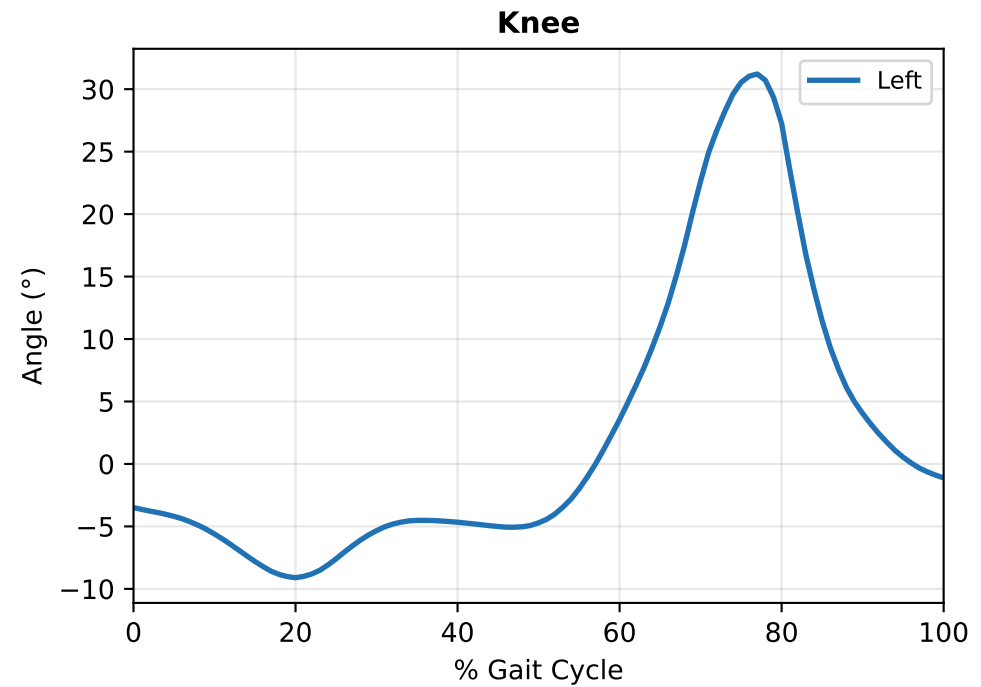


Bilateral Comparison

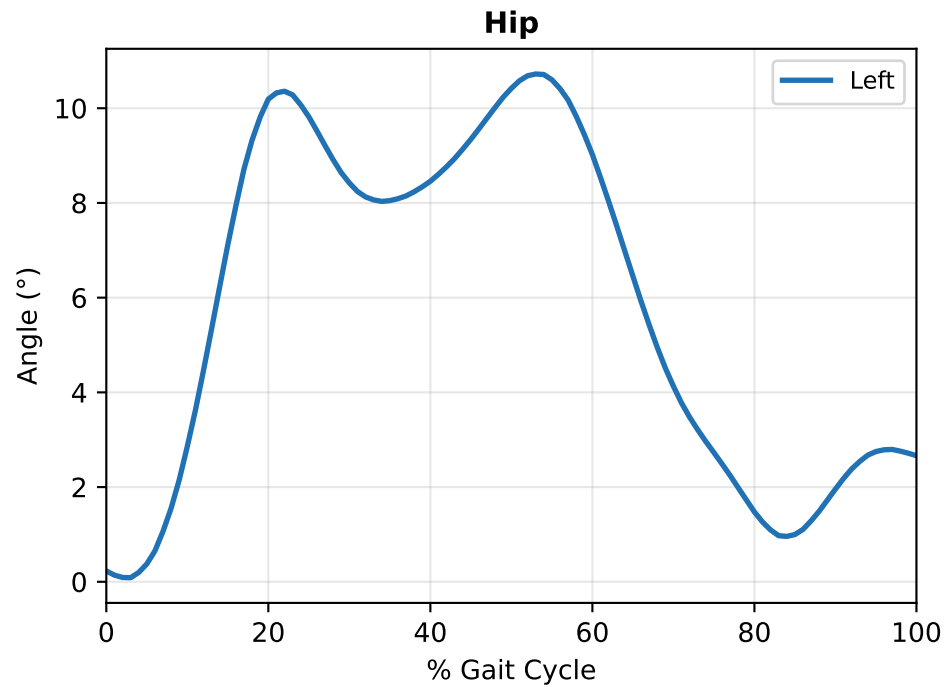


Source: 15-JUL-2026_pilot_ST_coronal_forward_1.MP4
FPS: 59.9
Cycles L: 1 R: 0
Model: mediapipe
Correction: 0.8

Clinical Statistics

Symmetry Indices

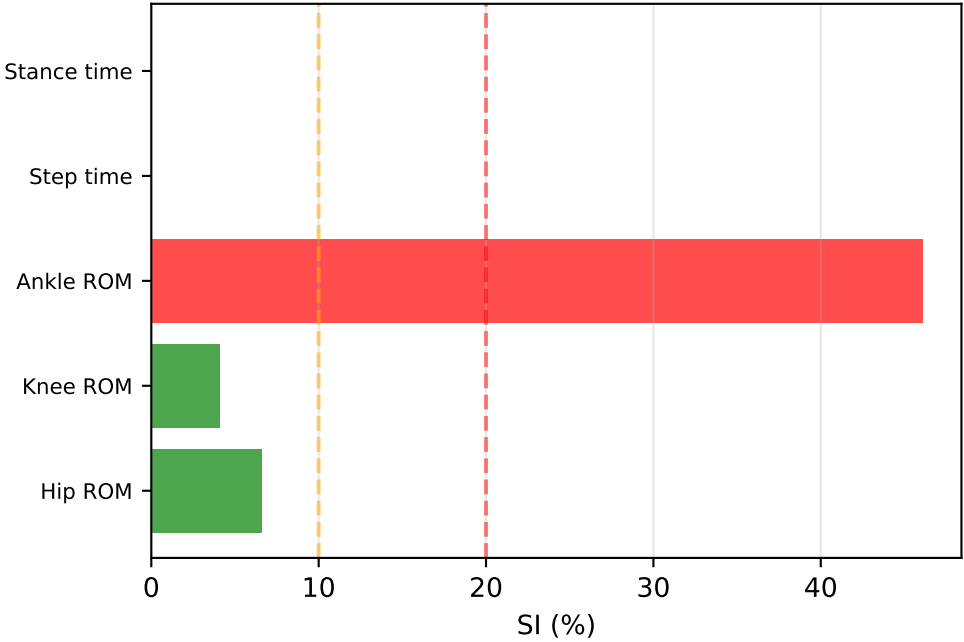

Variability (CV)

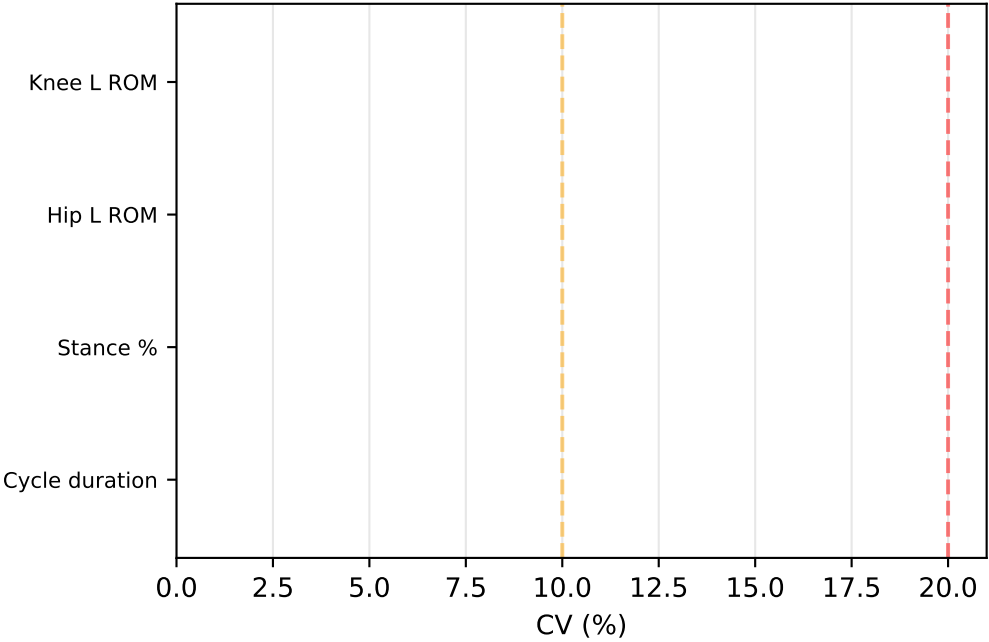

Temporal Parameters


Parameter	Value
Cadence	110.7 steps/min
Stride (mean)	1.084 s
Stride (SD)	0.0 s
Step (mean)	0.542 s
Stance L	53.8%
Stance R	None%
Double support	None%
Cycles	1

ALERTS


* Asymmetry ankle: SI=46.1% (>20%)

Trunk and Pelvis Analysis (Sagittal View)

Trunk Inclination — ROM: 3.2° Mean: 0.7° ± 0.8°

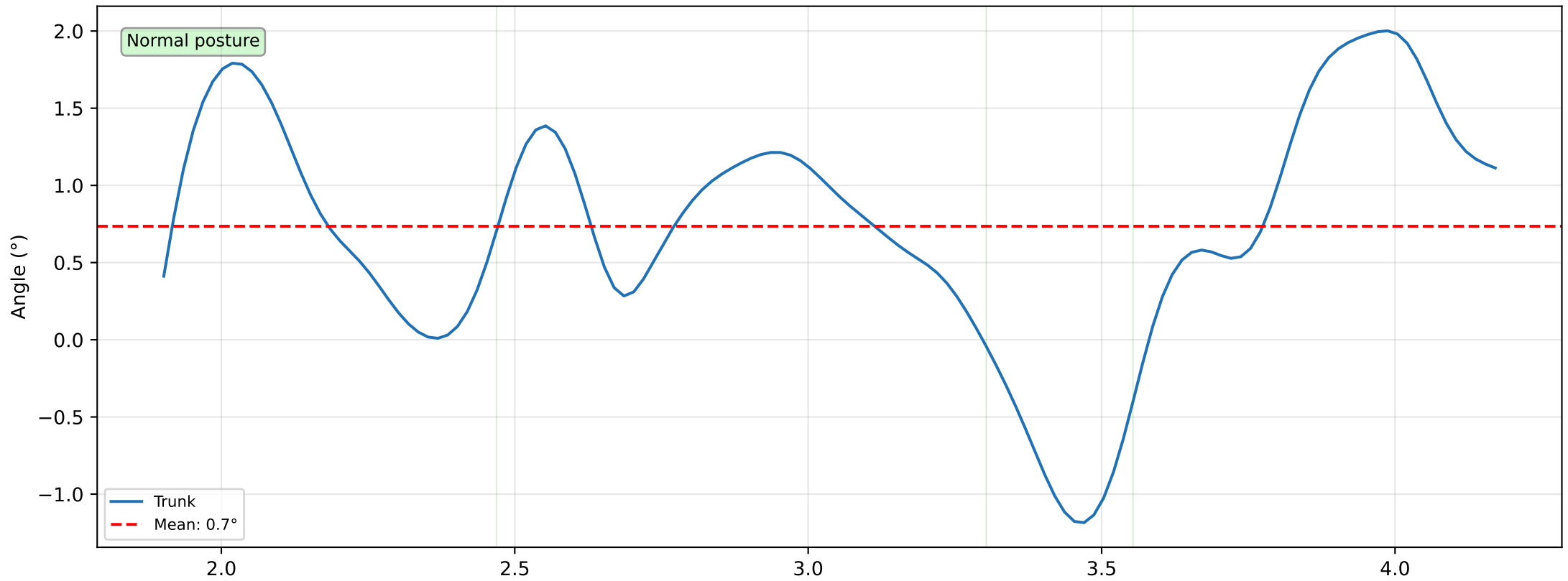

Pelvis Tilt — Requires Coronal View

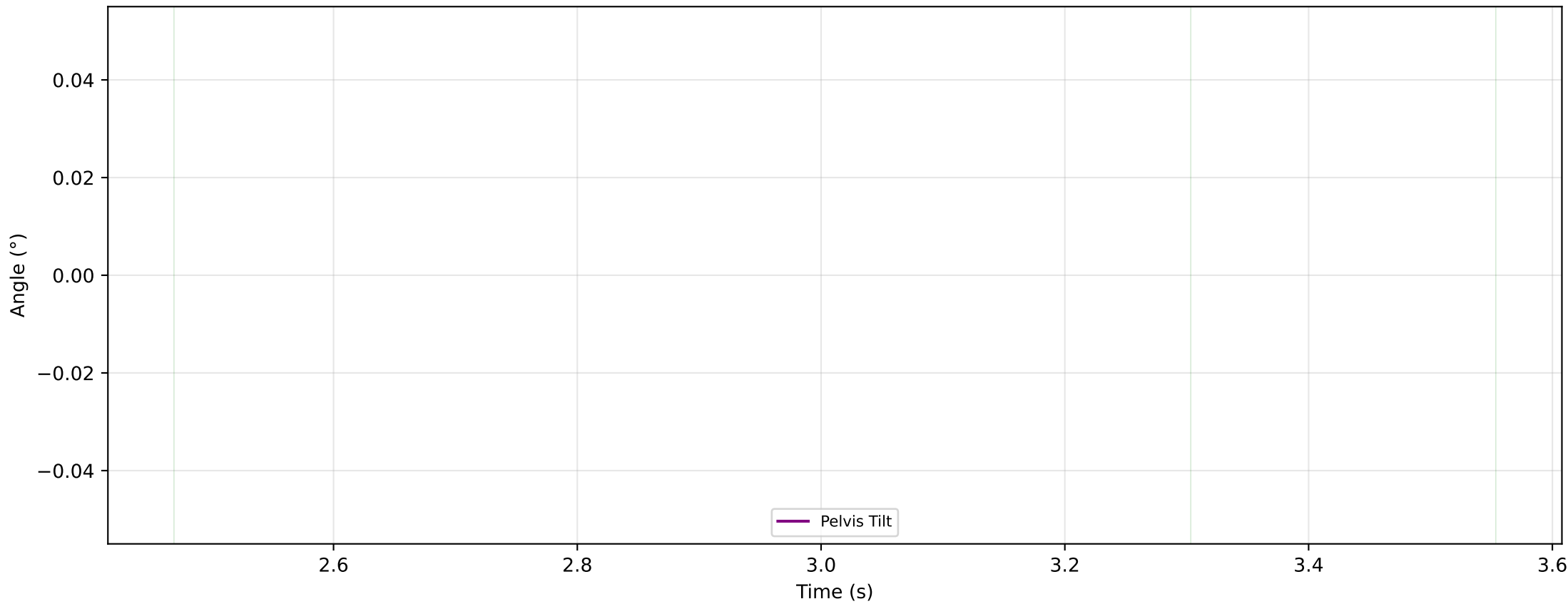


Pathology Detection


Pattern	Side	Severity	Value	Confidence	Description
Spastic	Left	mild	31.2	44%	Reduced swing knee flexion (31.2 deg, normal: 60-70)

Normalized Cycles — Left (n=1)


Hip — ROM mean: 10.6°

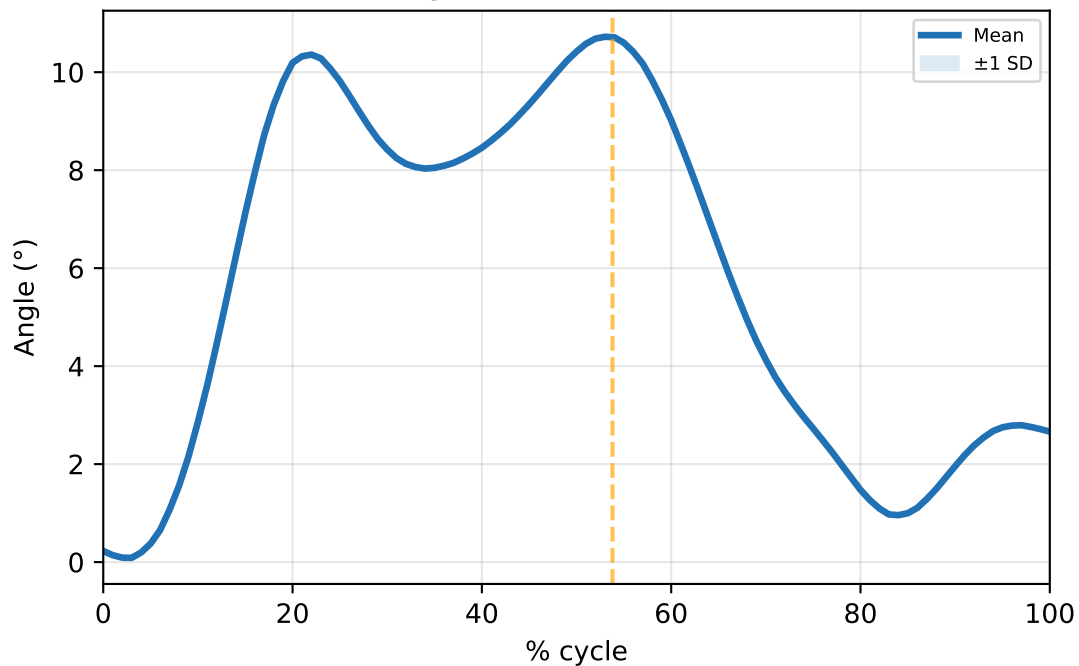


Knee — ROM mean: 40.3°

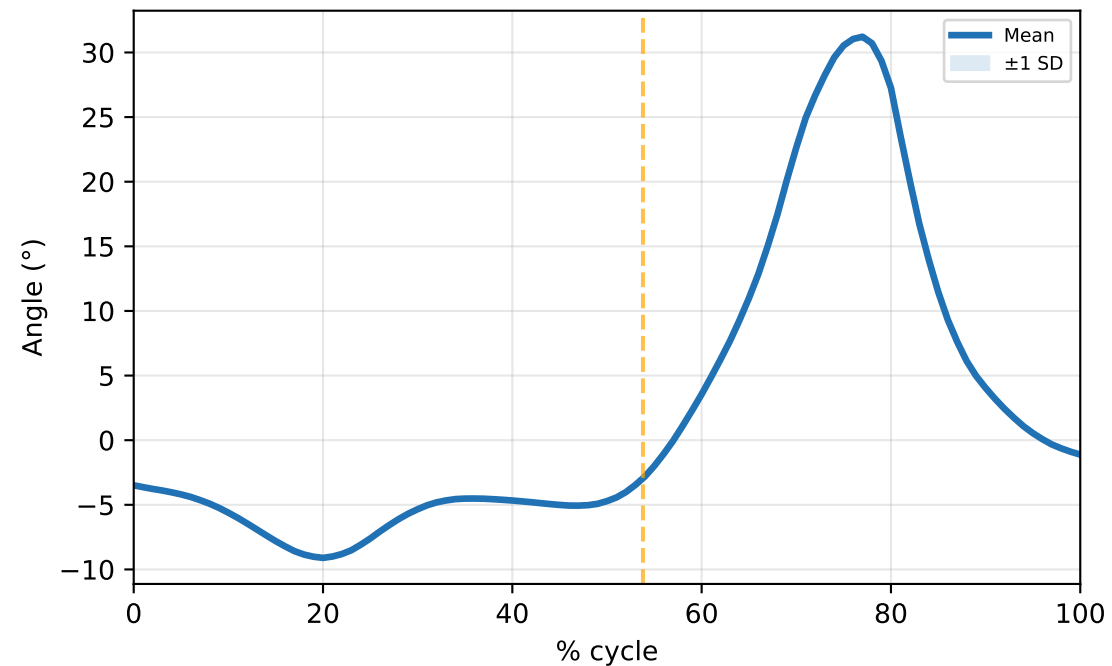

Ankle — ROM mean: 73.9°

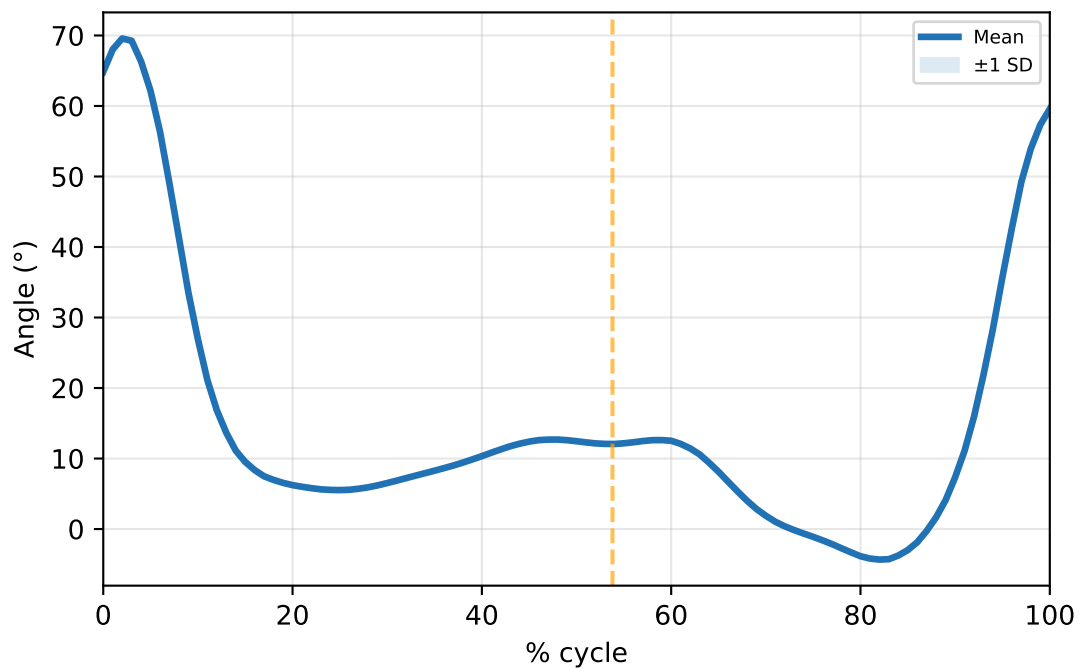


Trunk — ROM mean: 2.6°

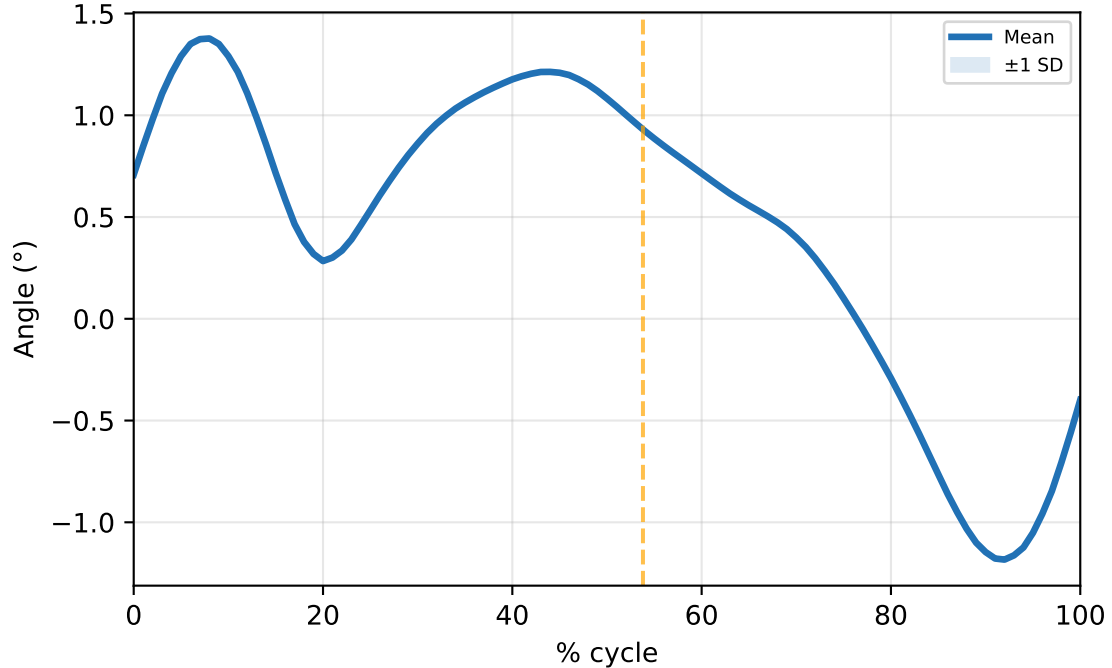


Normalized Cycles — Right (n=0)


Hip

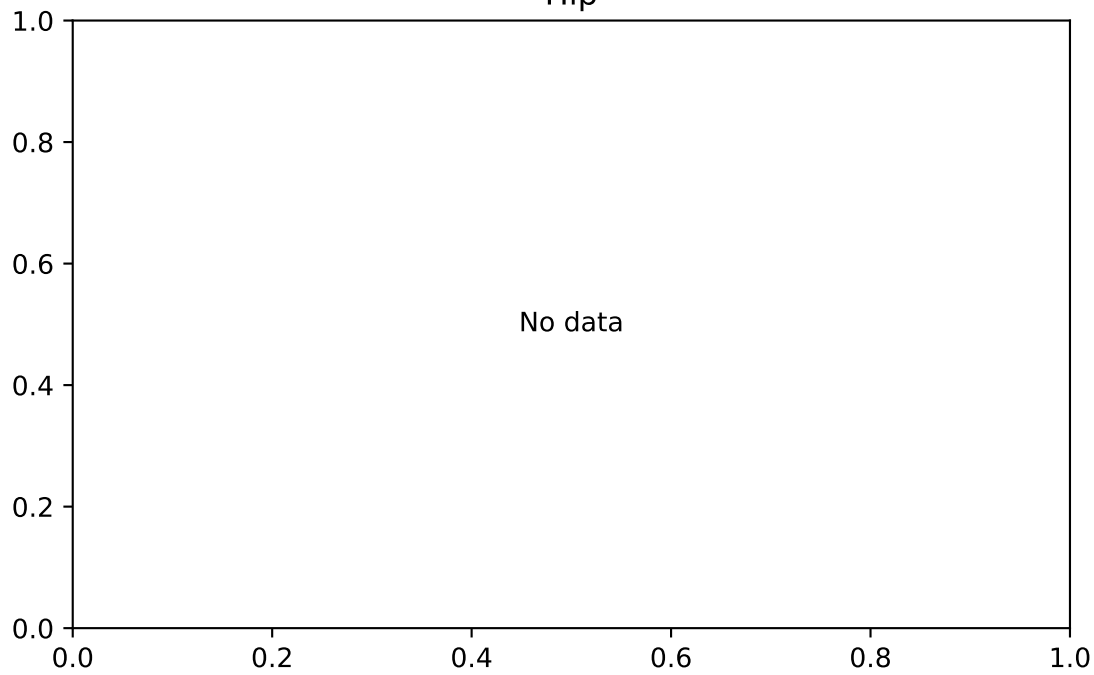


Knee

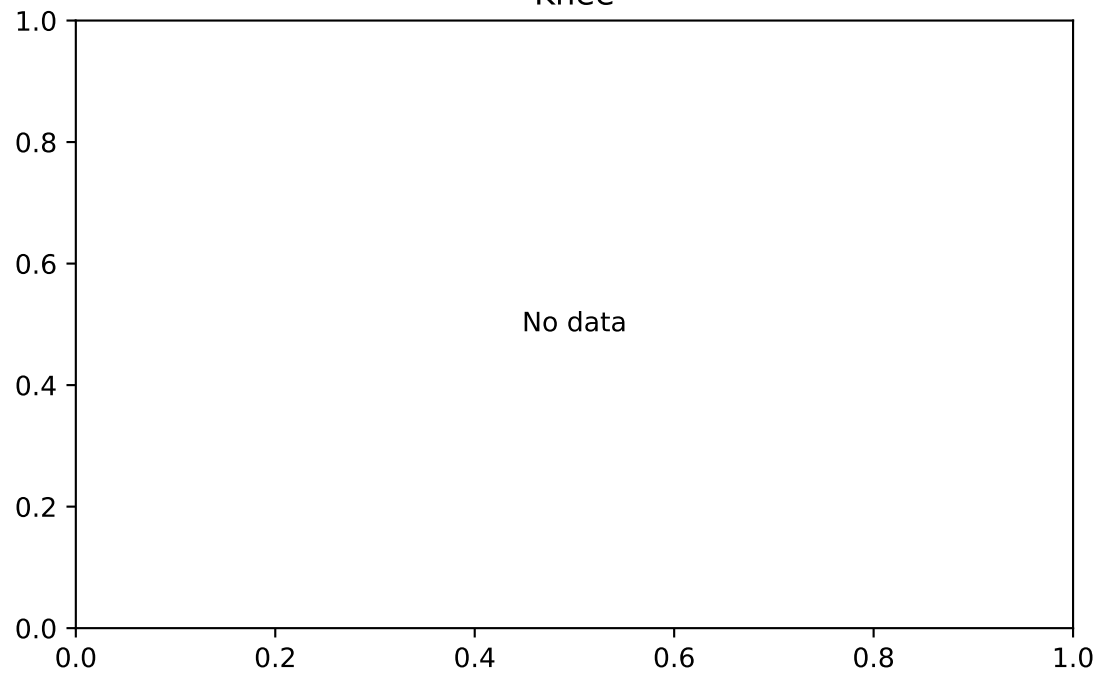

Ankle

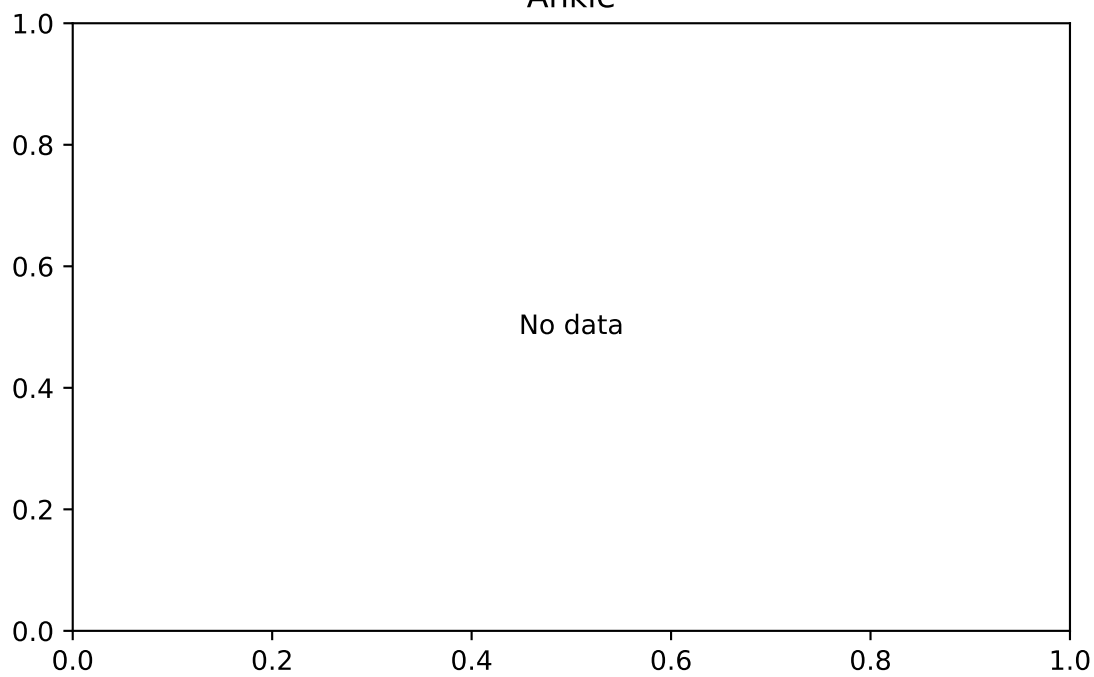


Trunk

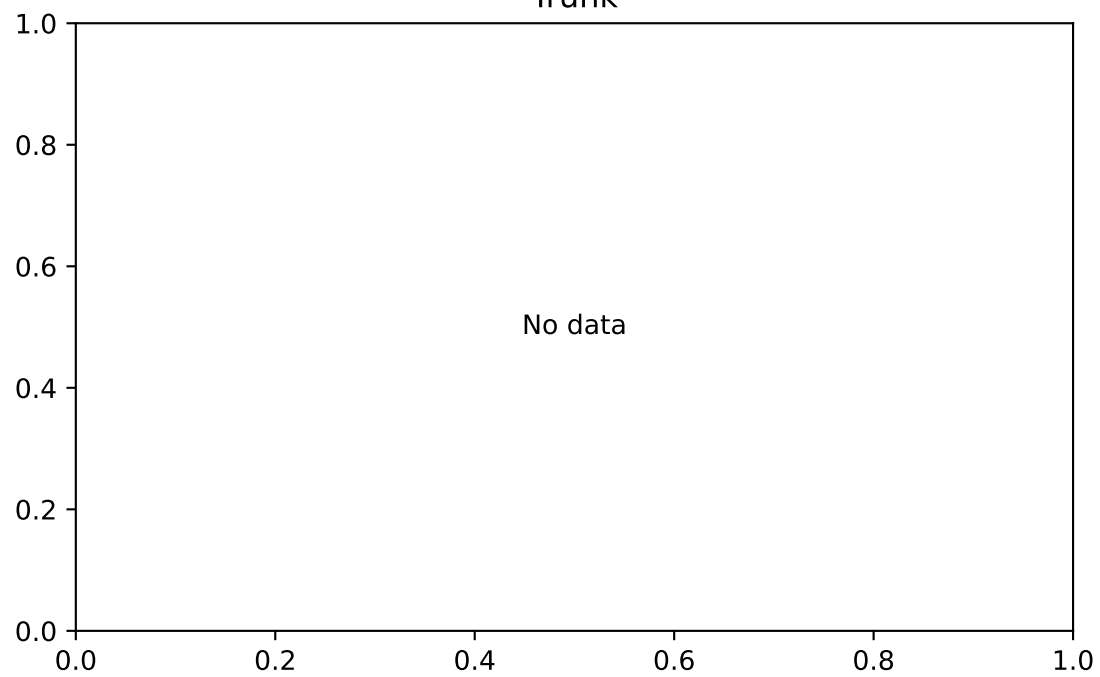


Normative Comparison


Hip — vs Normative (adult)

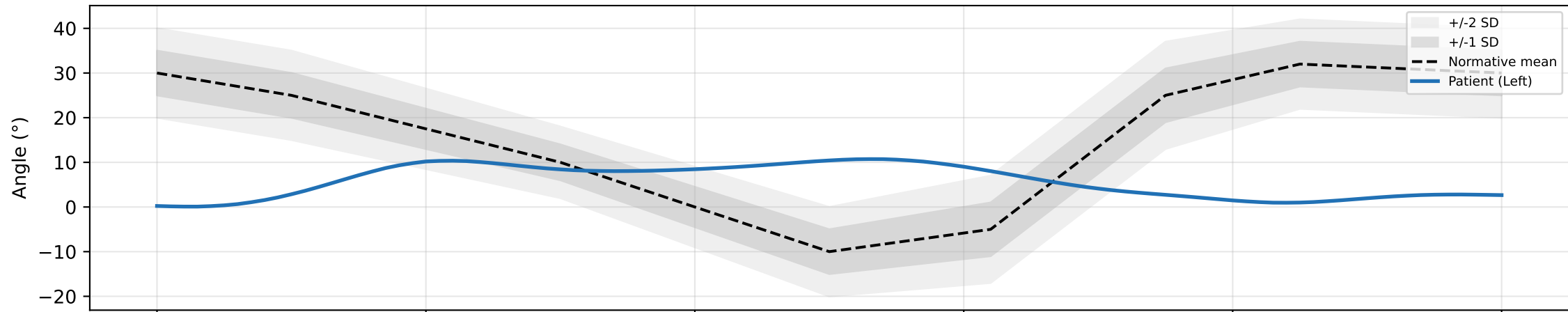


Knee — vs Normative (adult)

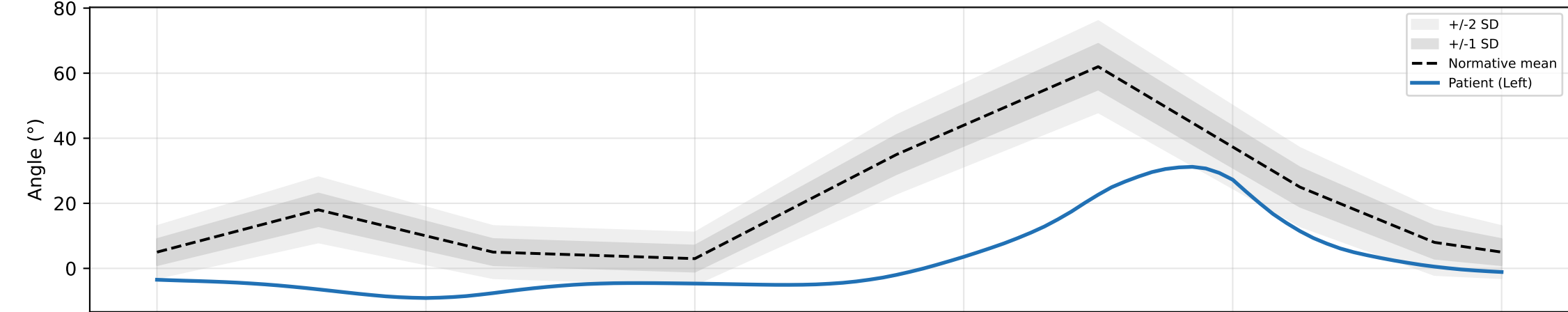

Ankle — vs Normative (adult)

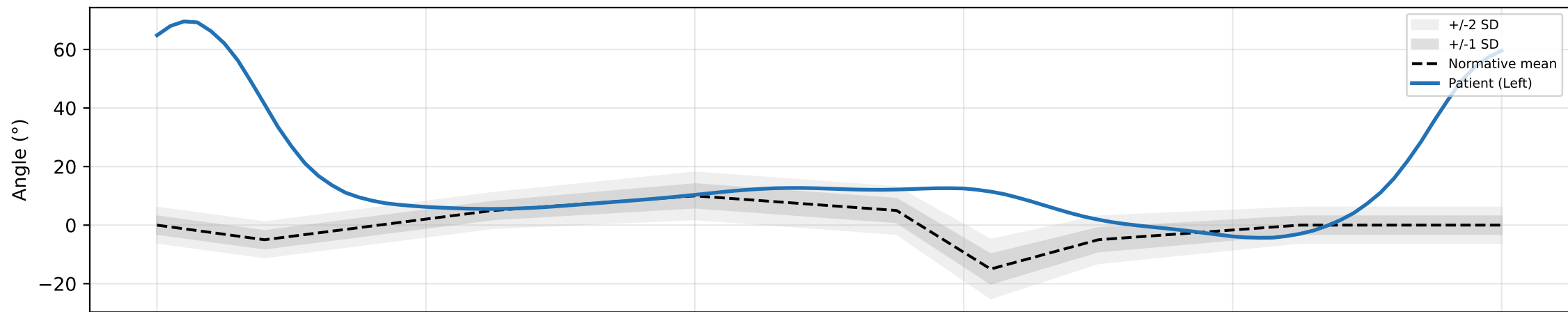


Trunk — vs Normative (adult)

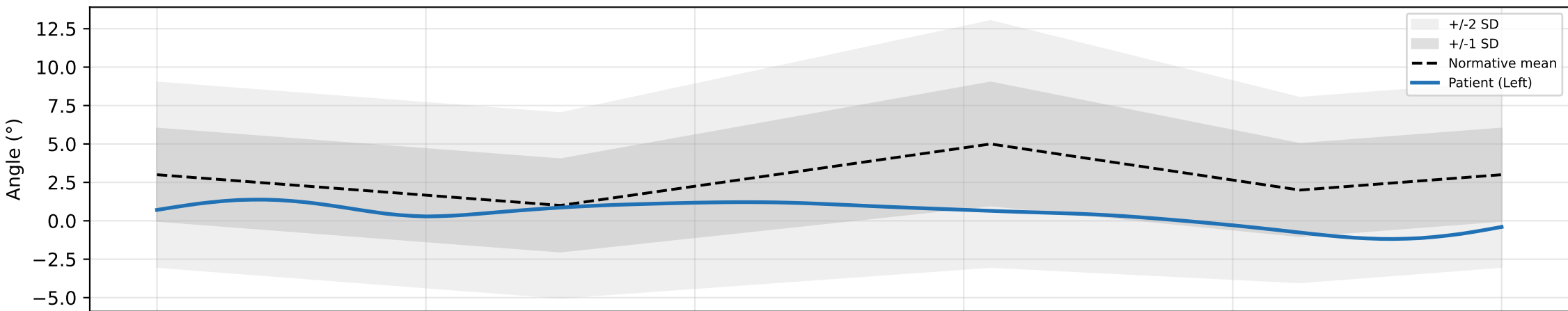


% Gait Cycle

GVS / MAP Profile

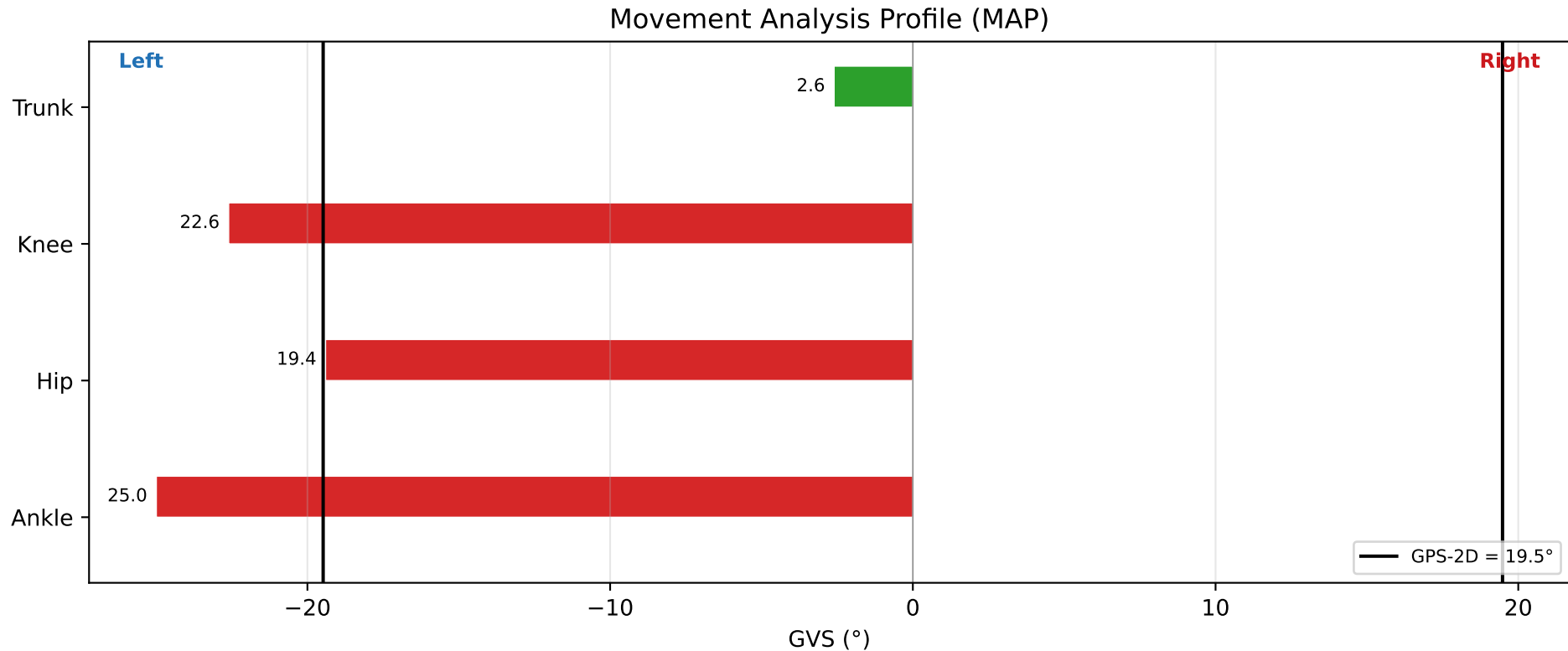


Quality Dashboard


Detection Rate per Landmark

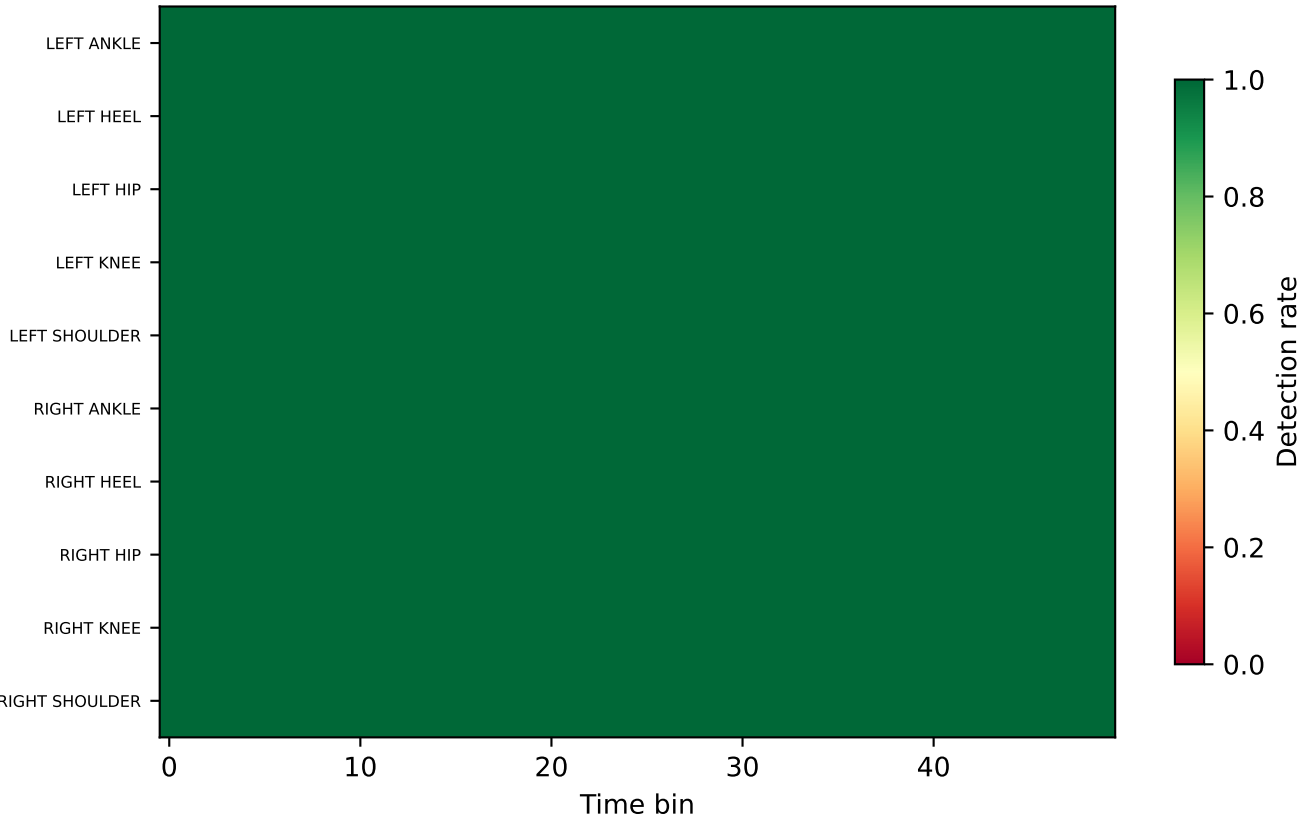

Per-frame Confidence

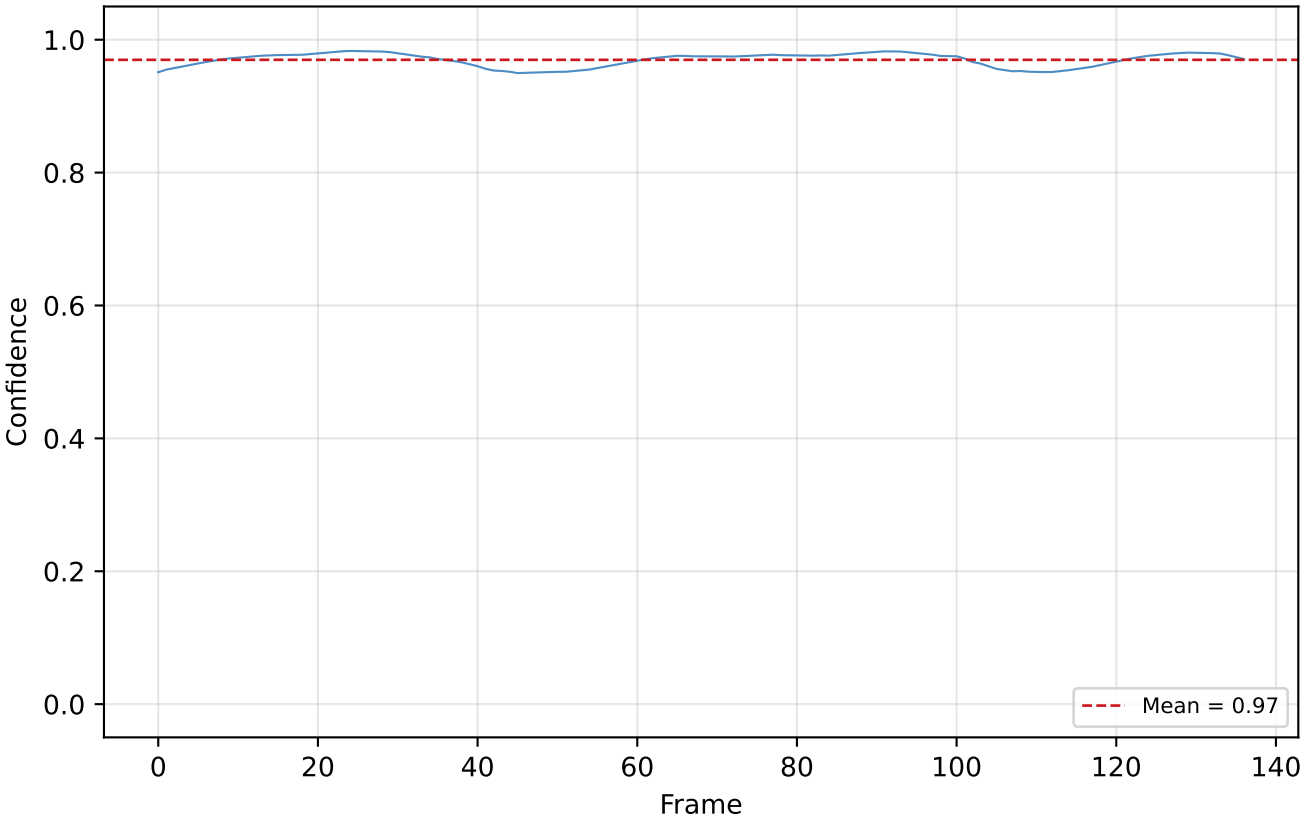

Gap Locations (0 gap frames)

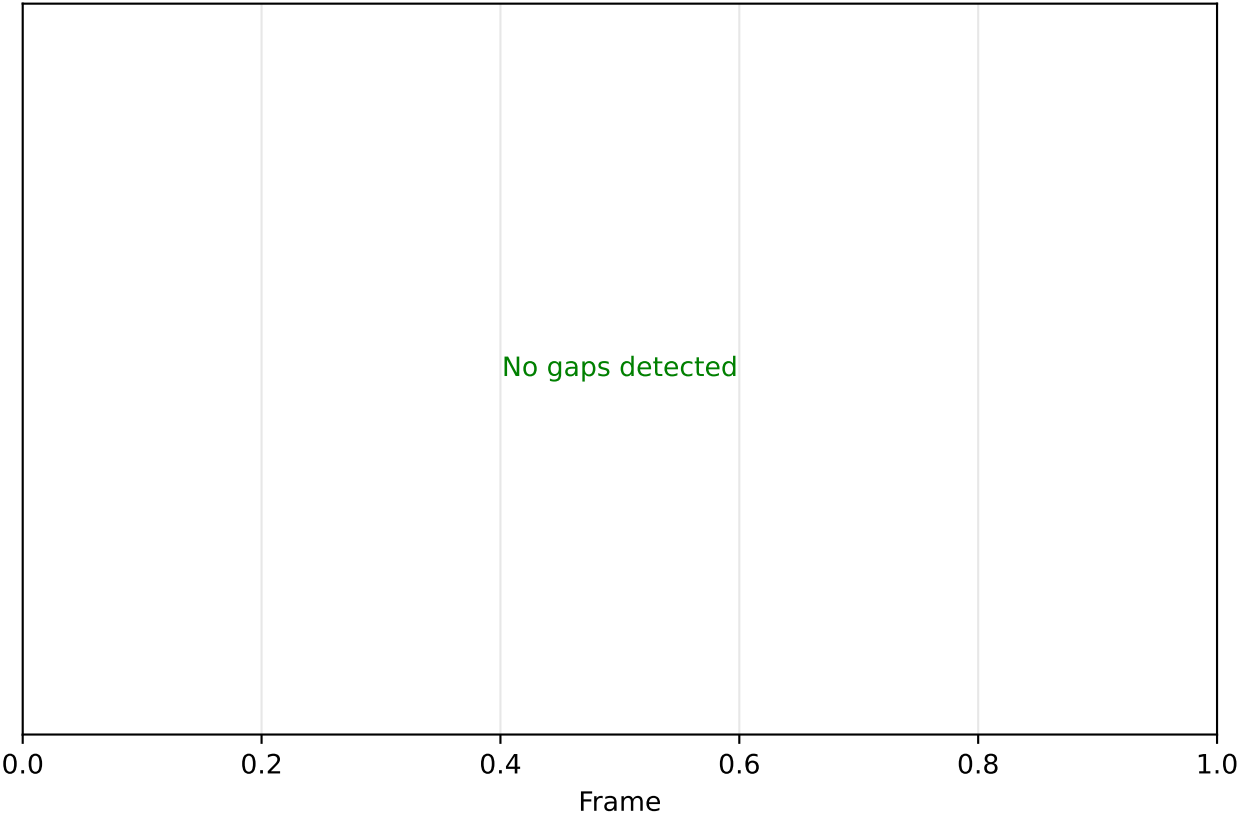

Overall Quality Score

99

Quality: GOOD

Detection rate: 100.0%

Mean confidence: 0.97

Gap percentage: 0.0%

Jitter score: 0.0008

Detailed Report

FULL REPORT – GAIT ANALYSIS

Source: /Users/adanmoallemi/Documents/MyoGait/videos/15-JUL-2026_pilot_ST_coronal_forward_1.MP4
FPS: 59.94005994005994 | Resolution: 1920x1080
Model: mediapipe | Correction: 0.8

SPATIO-TEMPORAL PARAMETERS

Cadence: 110.7 steps/min (normal: 100-120)
Stride (mean): 1.084 +/- 0.0 s
Step (mean): 0.542 +/- 0.292 s
Stance L: 53.8% Stance R: None% (normal: ~60%)
Double support: None% (normal: ~20%)

JOINT RANGES OF MOTION (ROM)

LEFT:
Hip 0.1 a 10.7 (ROM: 10.6)
Knee -9.1 a 31.2 (ROM: 40.3)
Ankle -4.3 a 69.6 (ROM: 73.9)

SYMMETRY

hip_rom_si 6.6%
knee_rom_si 4.1%
ankle_rom_si 46.1%
step_time_si N/A%
overall_si 18.9%
(SI < 10% = excellent, 10-20% = good, >20% = asymmetry)

VARIABILITY

Cycle duration CV: 0.0%
Stance % CV: 0.0%

ALERTS

* Asymmetry ankle: SI=46.1% (>20%)

Normal reference values:

Hip ROM: 40-50° | Knee ROM: 60-70° | Ankle ROM: 30-40°

Generated by myogait v0.6.0