

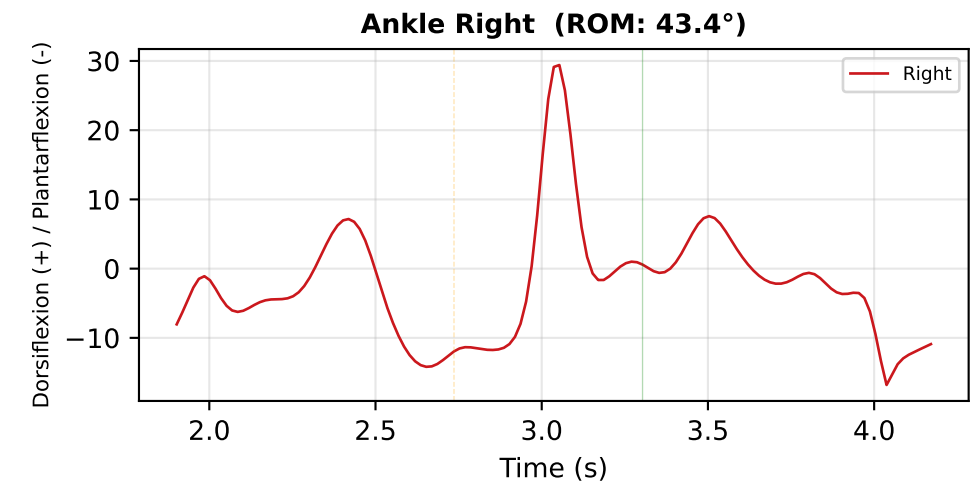
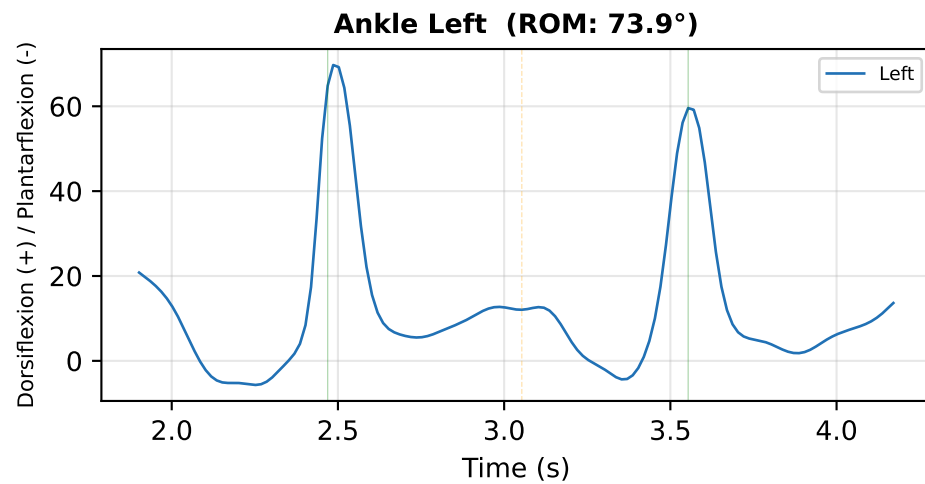
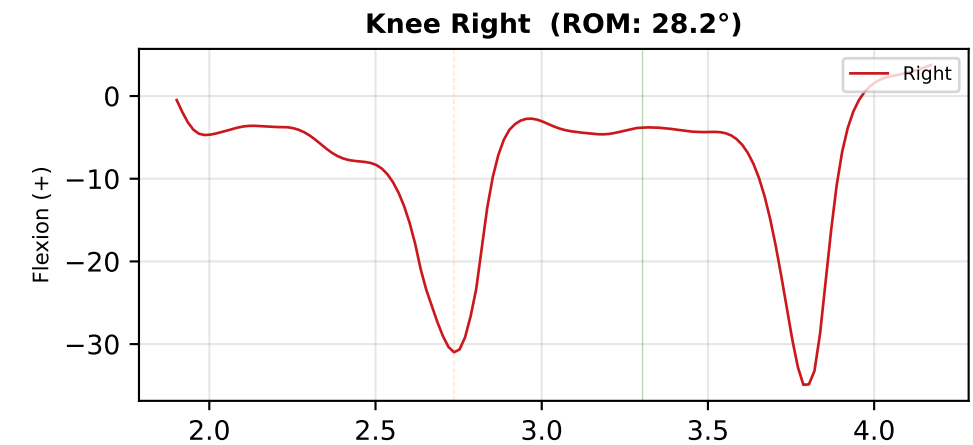
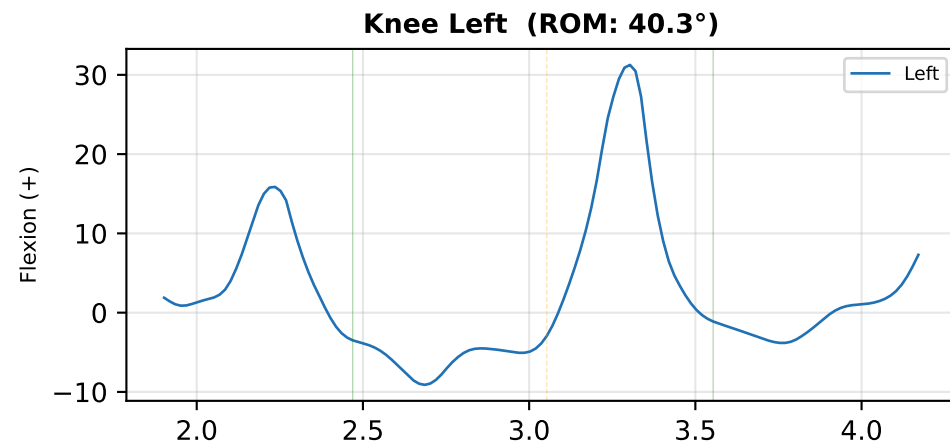
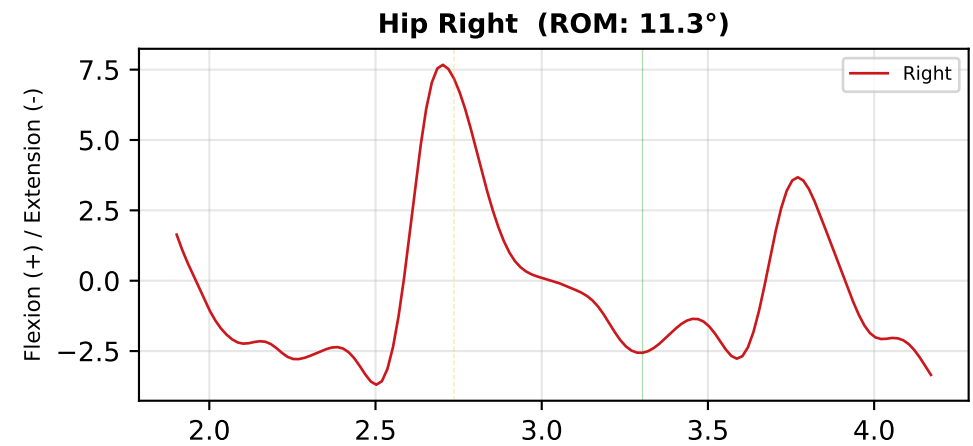
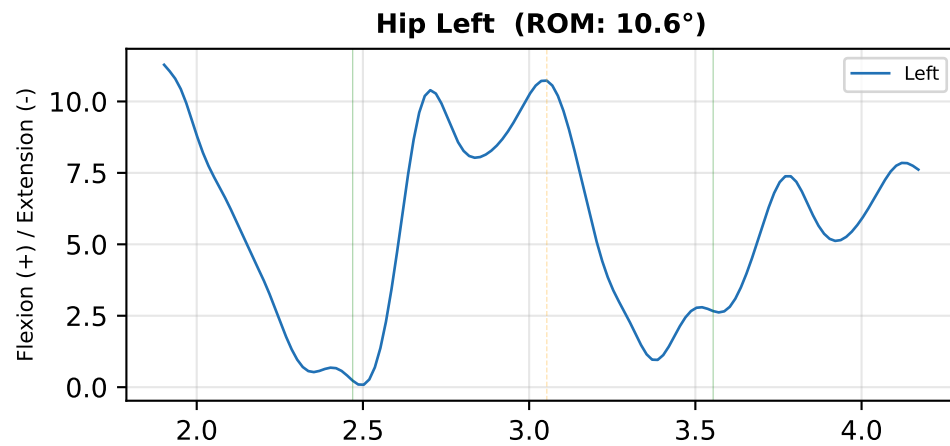
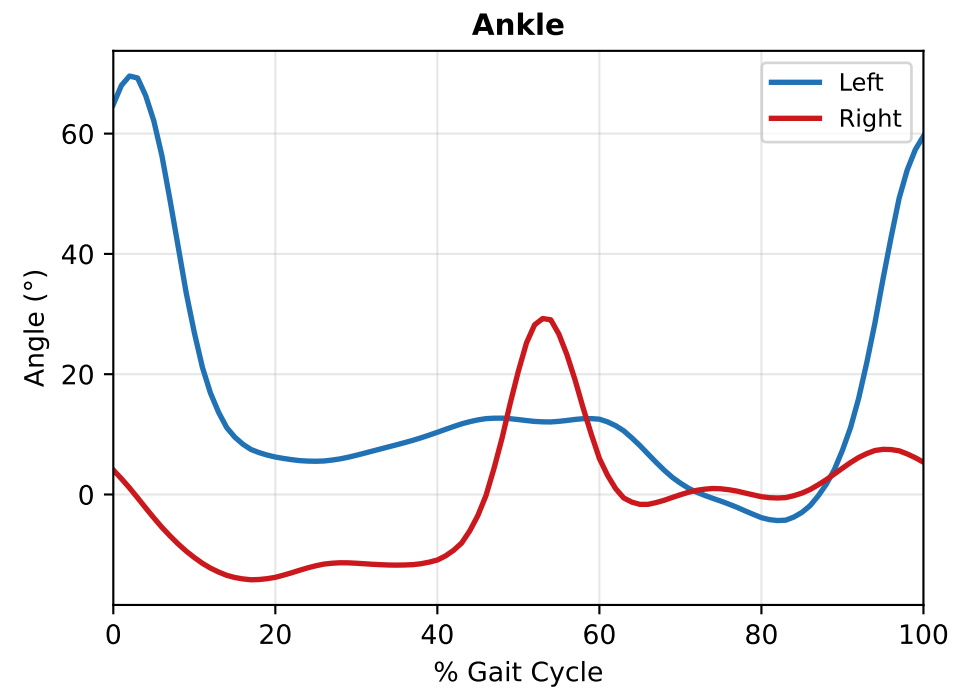


Overview — Joint Angles


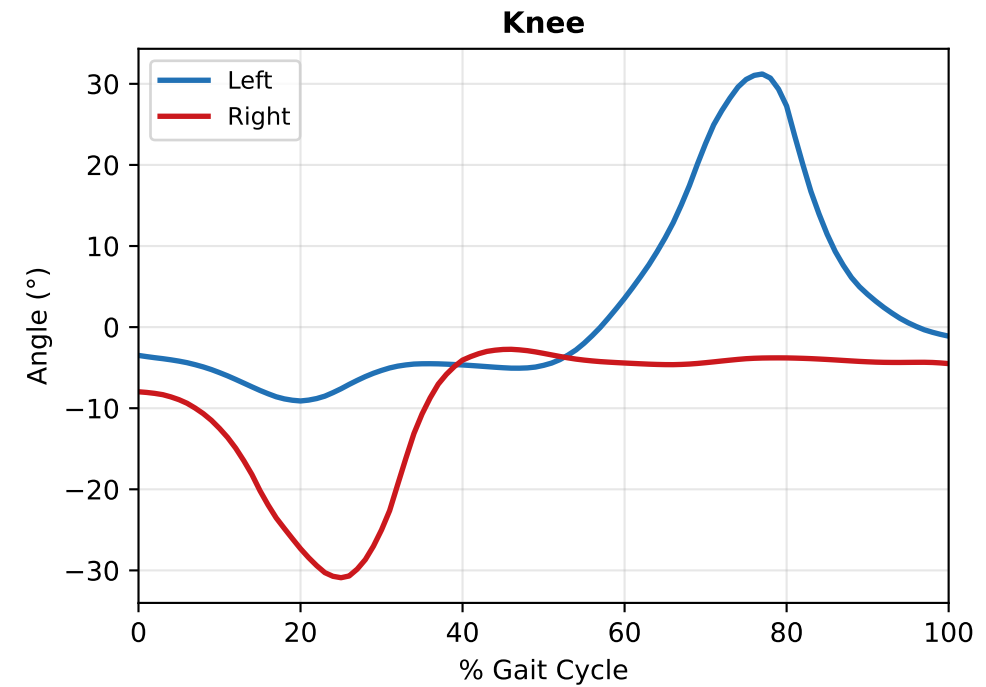


Bilateral Comparison

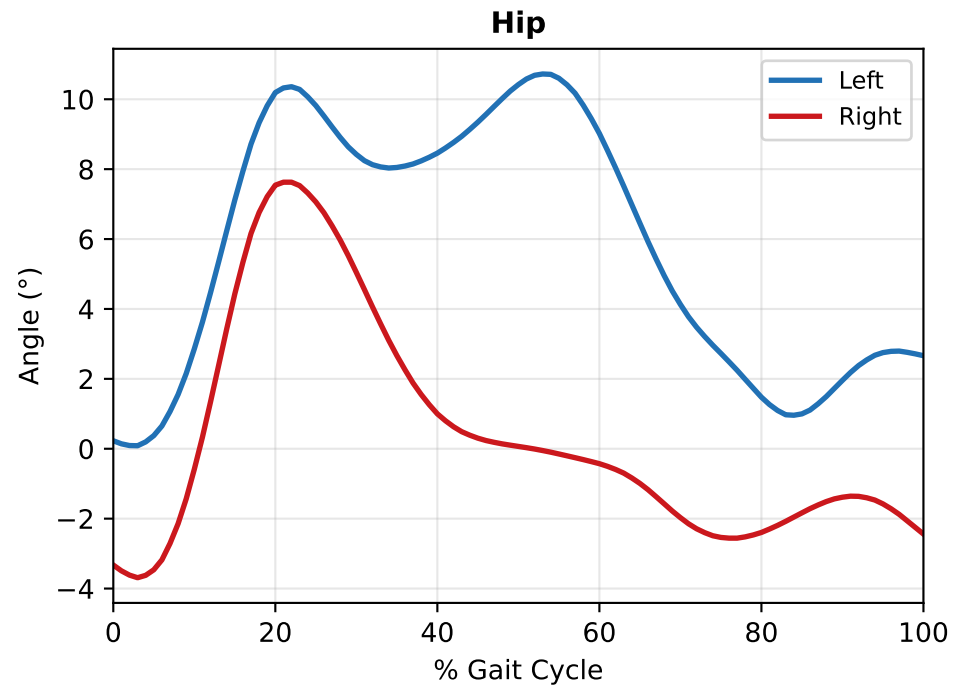

Source: 15-JUL-2026_pilot_ST_coronal_forward_1.MP4
FPS: 59.9
Cycles L: 1 R: 1
Model: mediapipe
Correction: 0.8

Clinical Statistics

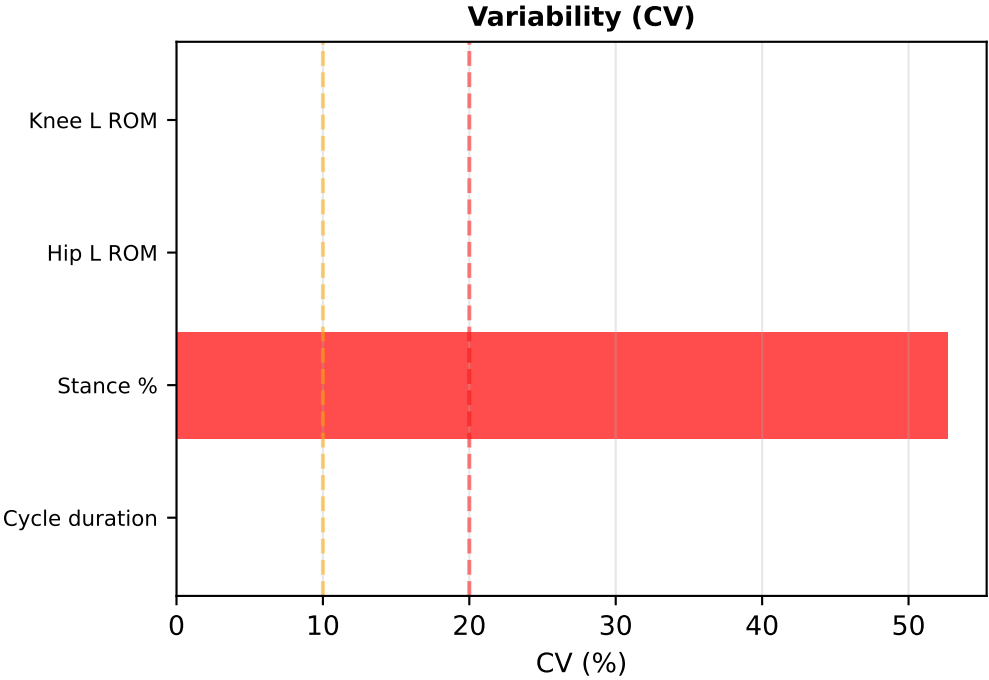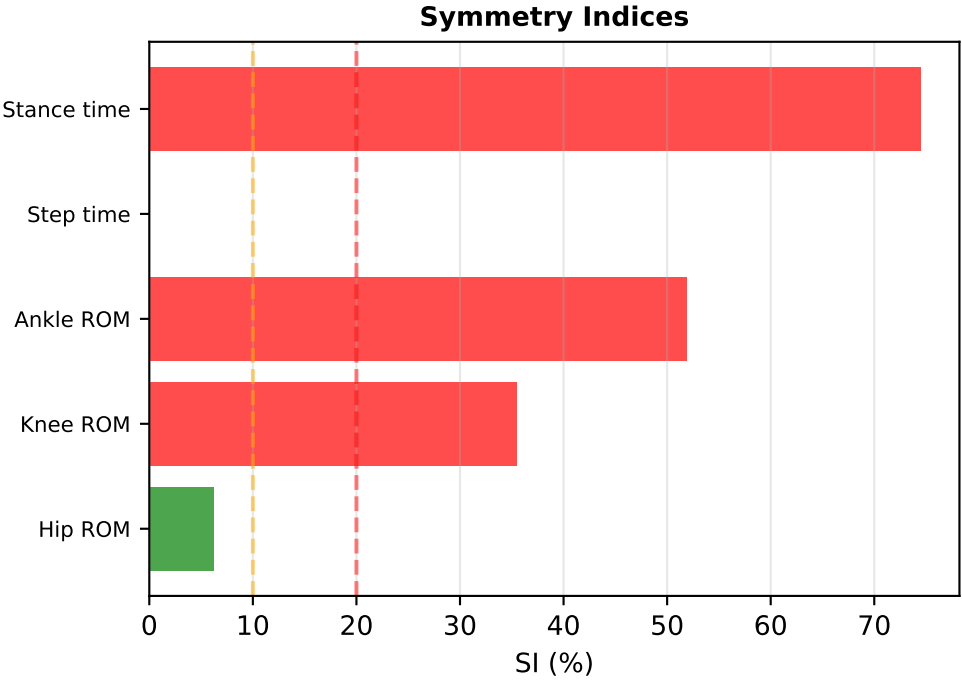

Temporal Parameters

Parameter	Value
Cadence	110.7 steps/min
Stride (mean)	1.084 s
Stride (SD)	0.0 s
Step (mean)	0.542 s
Stance L	53.8%
Stance R	None%
Double support	None%
Cycles	2

ALERTS


* Asymmetry knee: SI=35.5% (>20%)

* Asymmetry ankle: SI=51.9% (>20%)


* Asymmetry stance_time: SI=74.5% (>20%)

Trunk and Pelvis Analysis (Coronal View)

Trunk Inclination — ROM: 3.2° Mean: $0.7^{\circ} \pm 0.8^{\circ}$

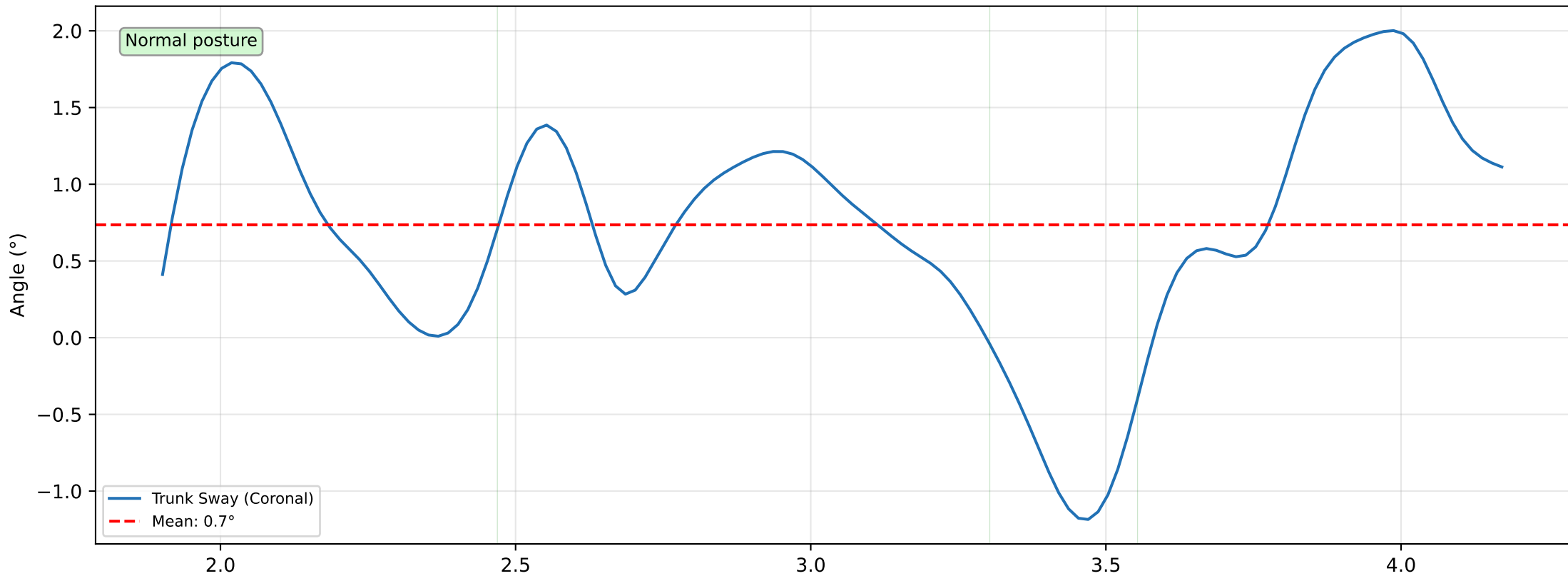

Pelvis List (Coronal) — ROM: 8.6° Mean: -0.1°

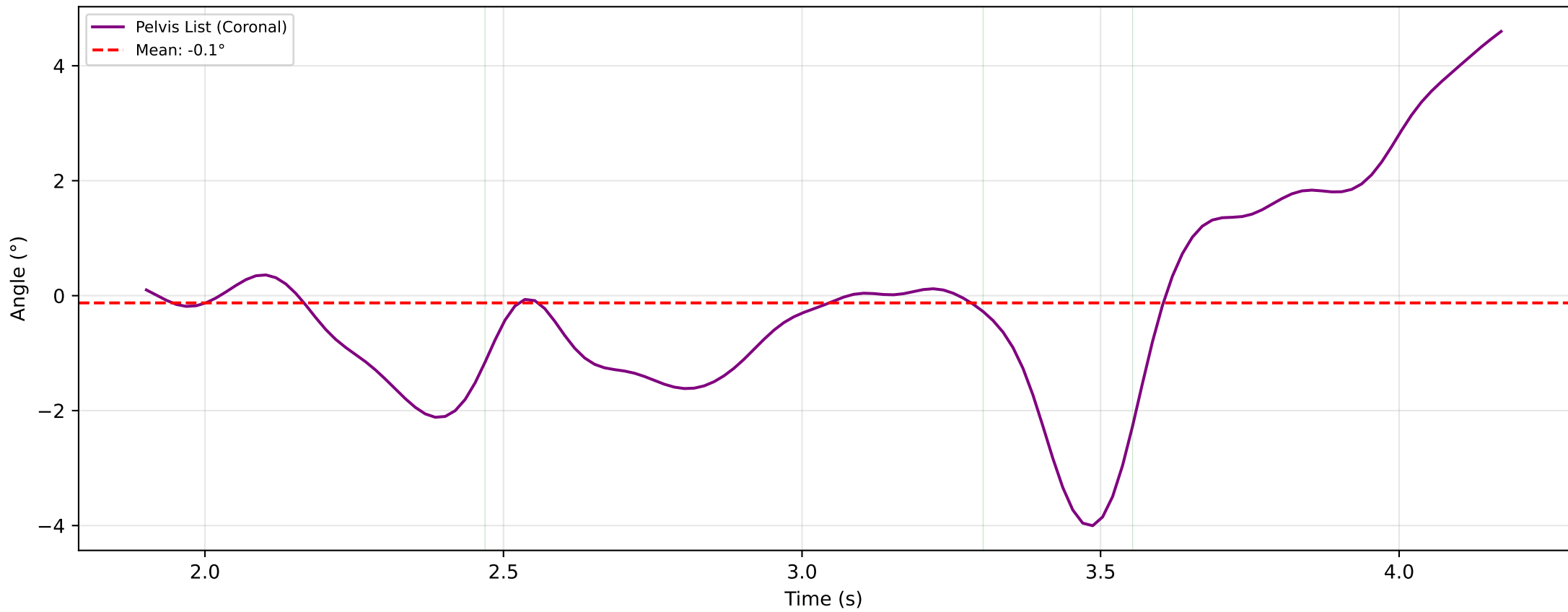


Pathology Detection


Pattern	Side	Severity	Value	Confidence	Description
Spastic	Left	mild	31.2	44%	Reduced swing knee flexion (31.2 deg, normal: 60-70)
Spastic	Right	moderate	-3.8	100%	Reduced swing knee flexion (-3.8 deg, normal: 60-70)

Normalized Cycles — Left (n=1)


Hip — ROM mean: 10.6°

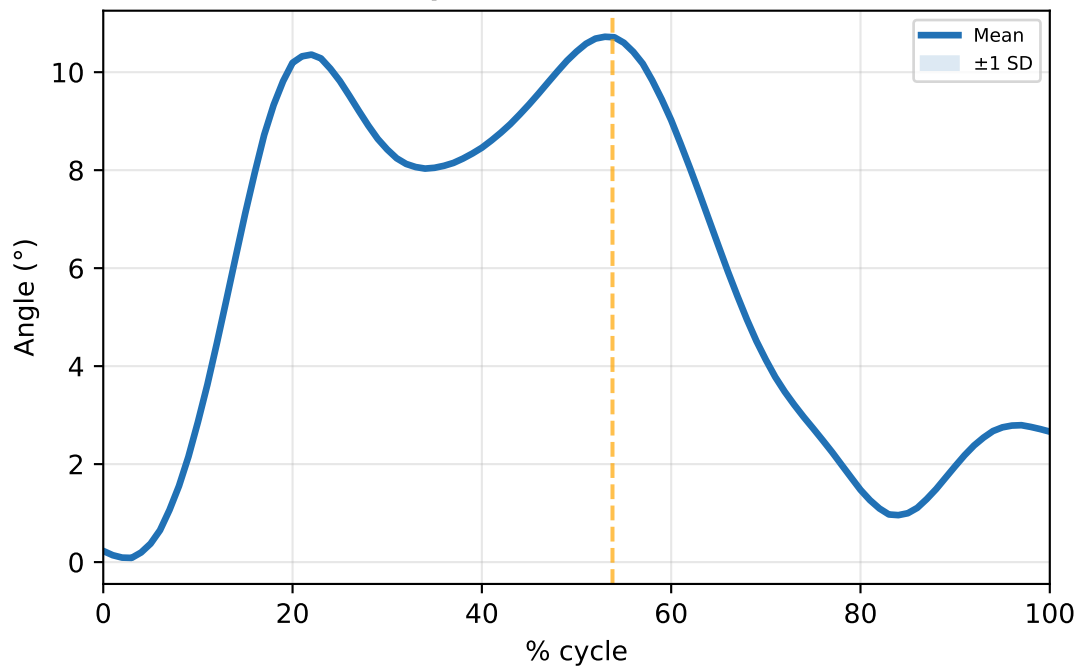


Knee — ROM mean: 40.3°

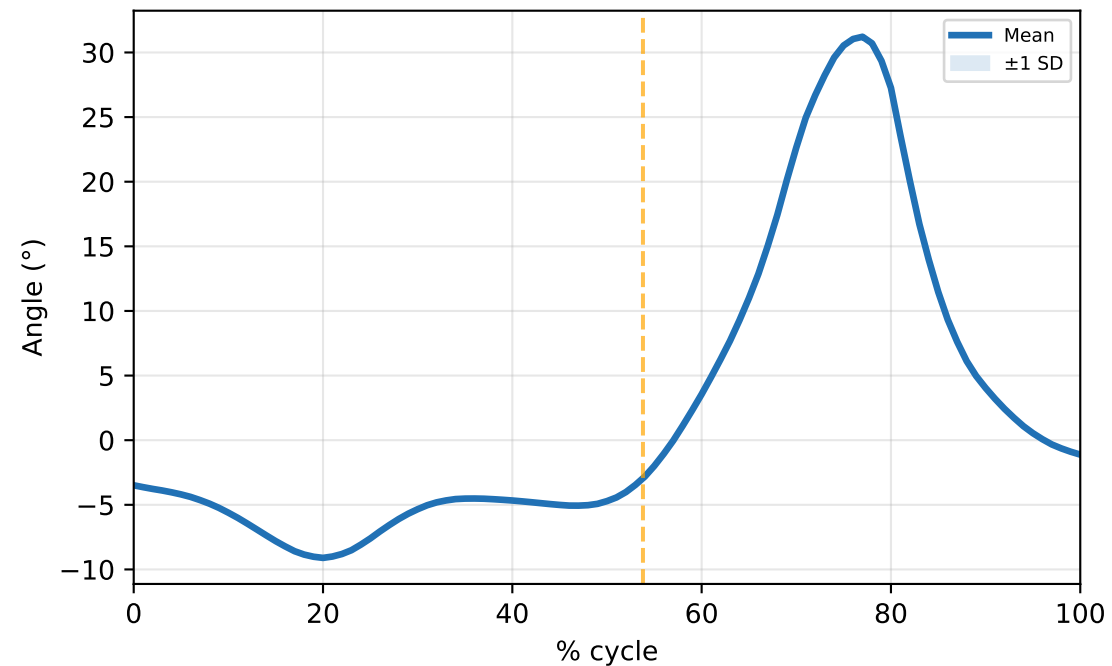

Ankle — ROM mean: 73.9°

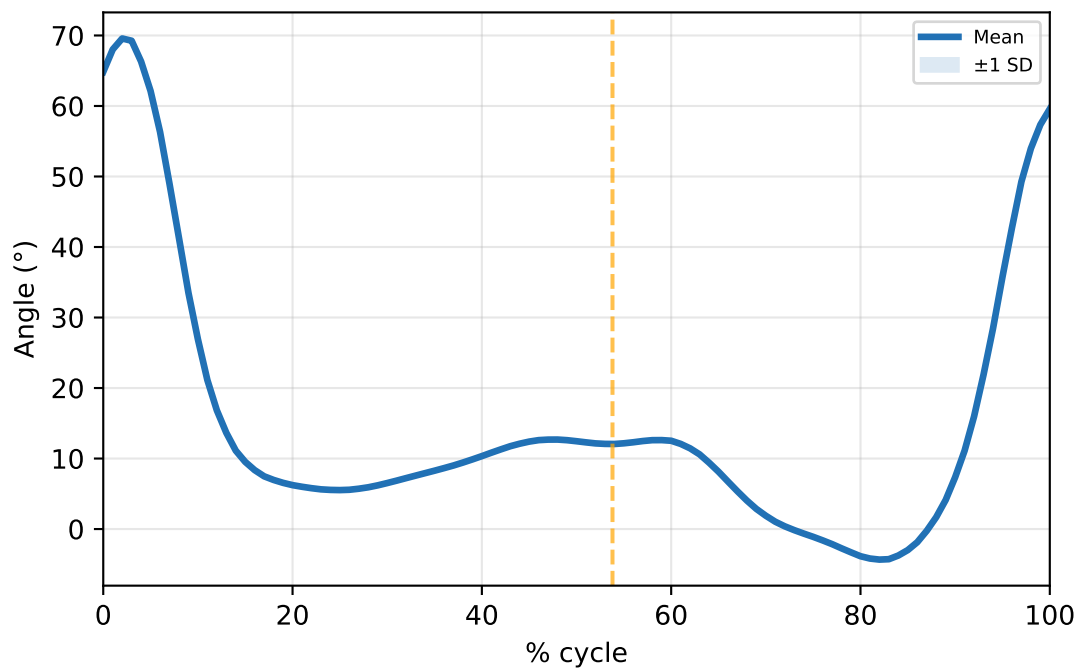


Trunk — ROM mean: 2.6°

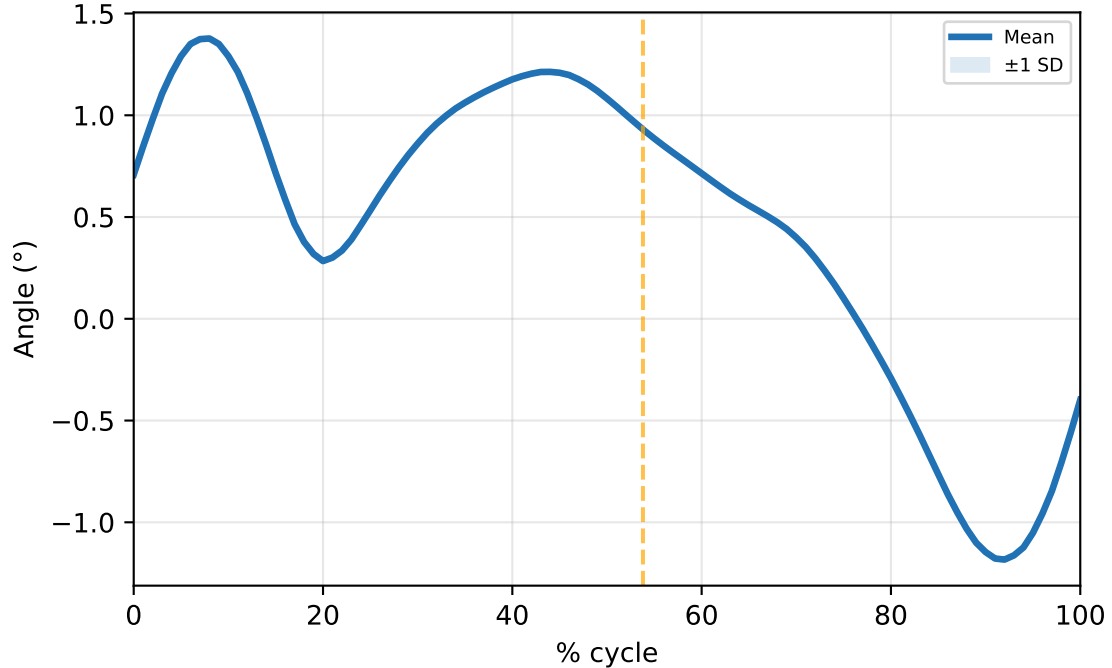


Normalized Cycles — Right (n=1)


Hip — ROM mean: 11.3°

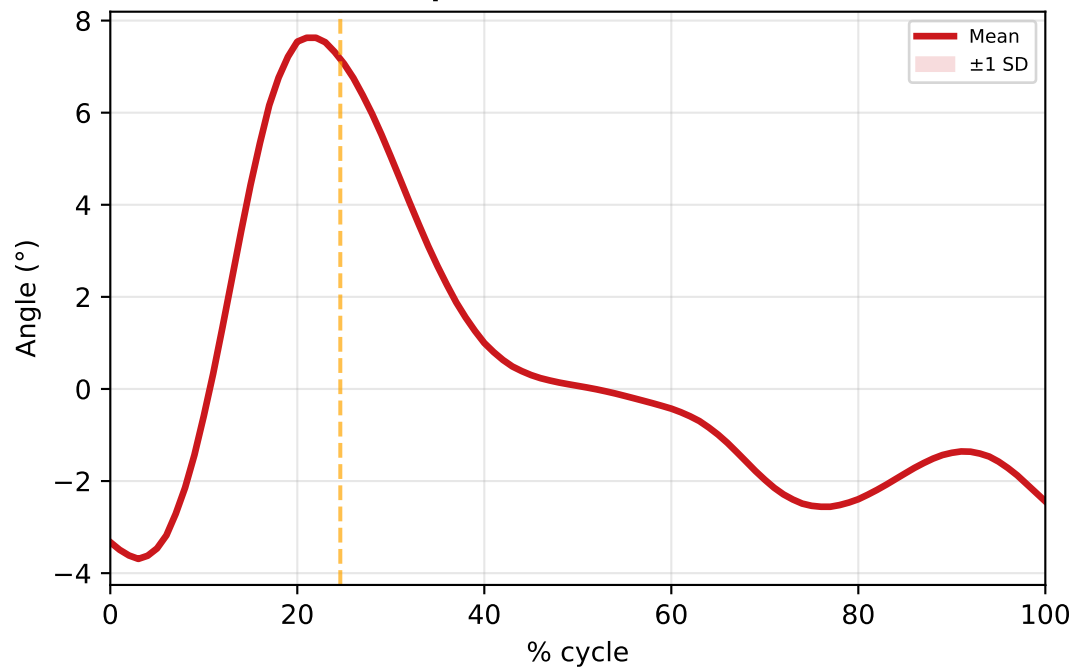


Knee — ROM mean: 28.2°

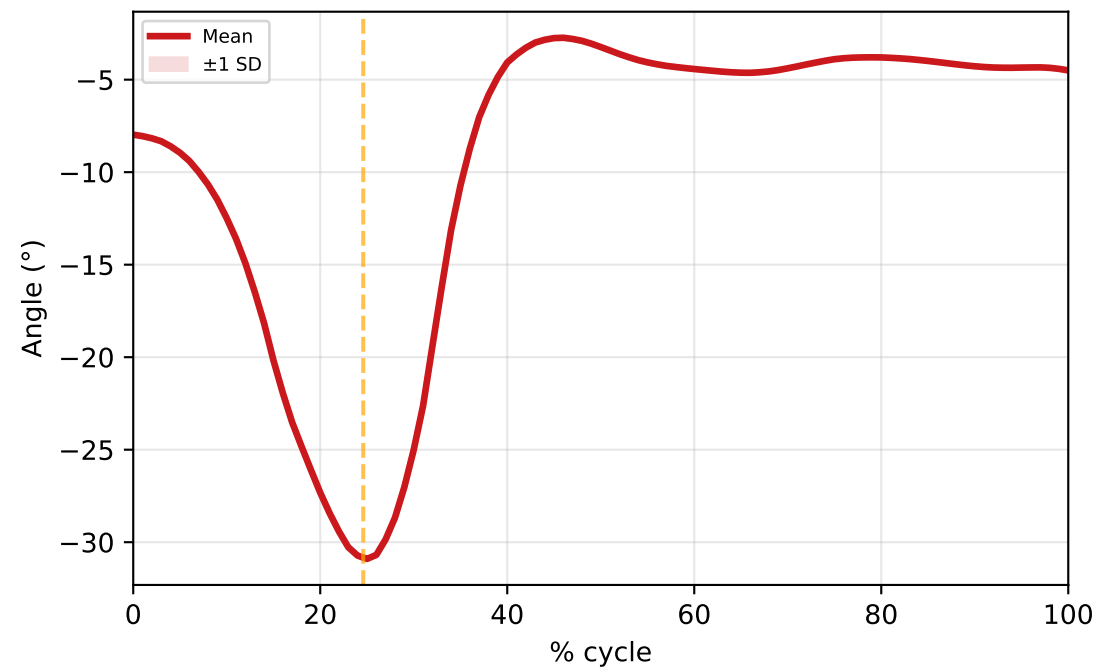

Ankle — ROM mean: 43.4°

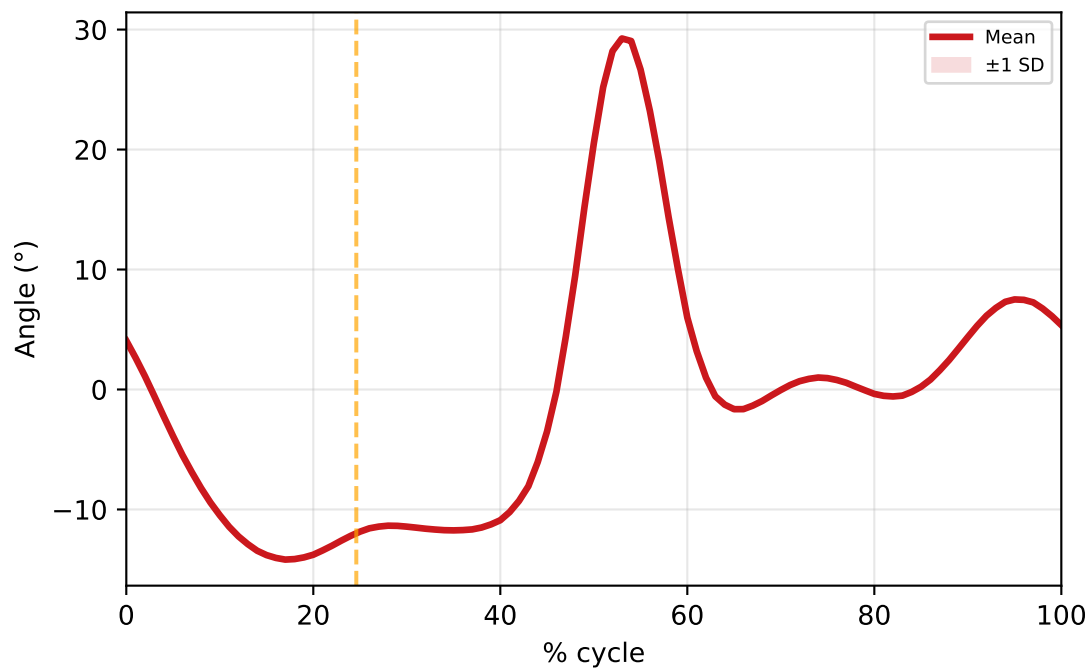


Trunk — ROM mean: 2.6°

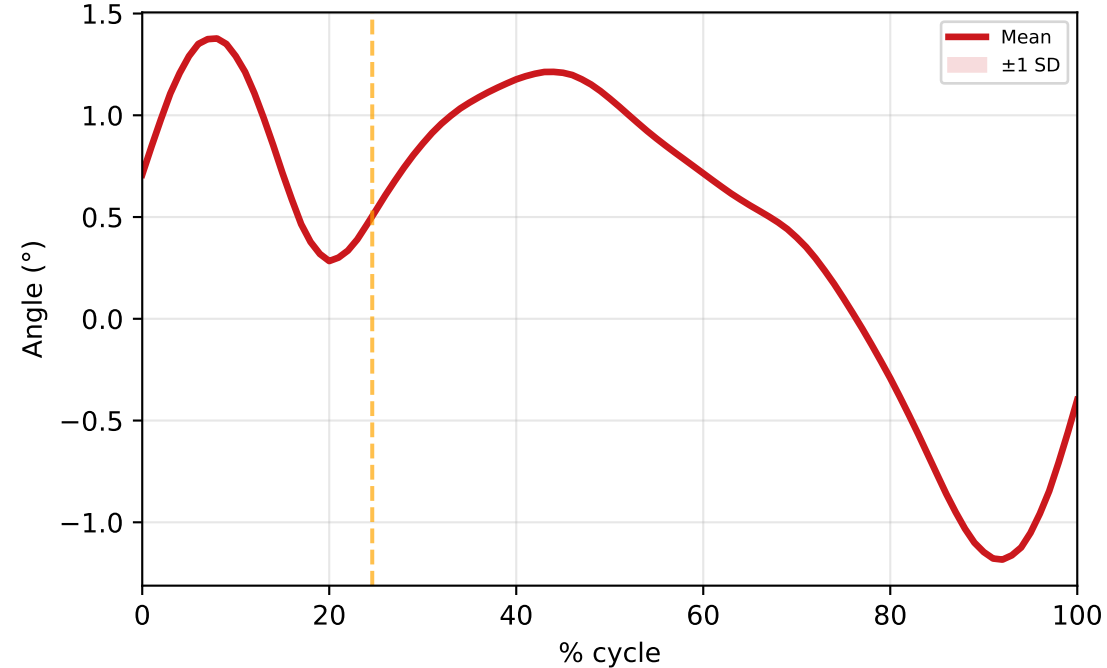


Normative Comparison


Hip — vs Normative (adult)

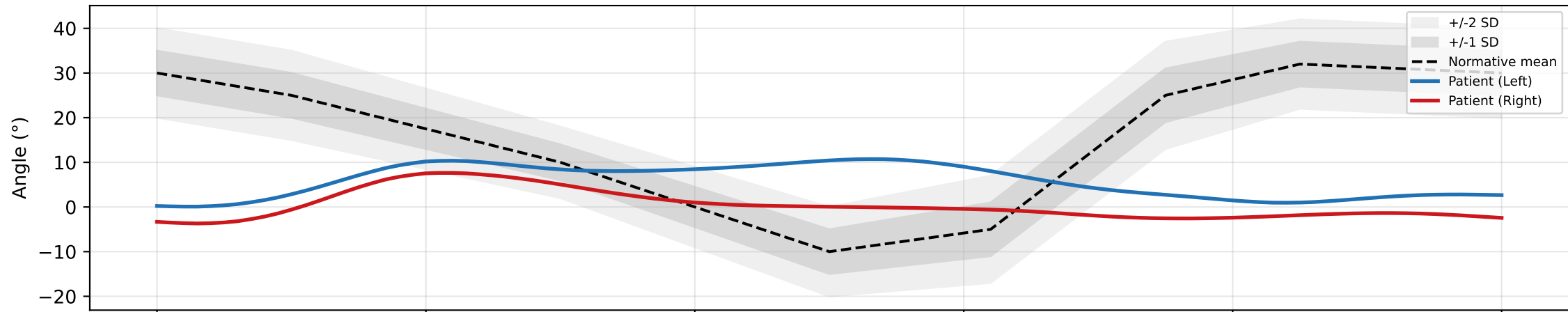


Knee — vs Normative (adult)

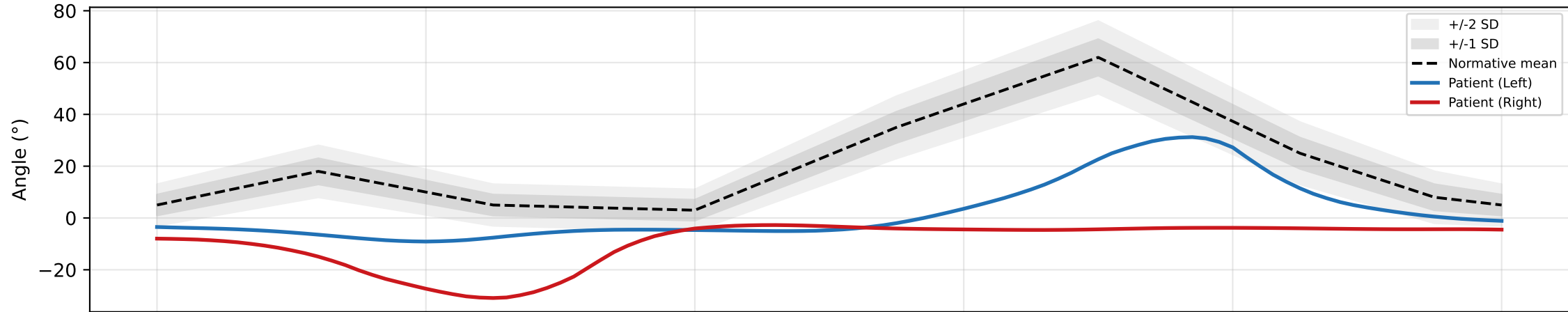


Ankle — vs Normative (adult)

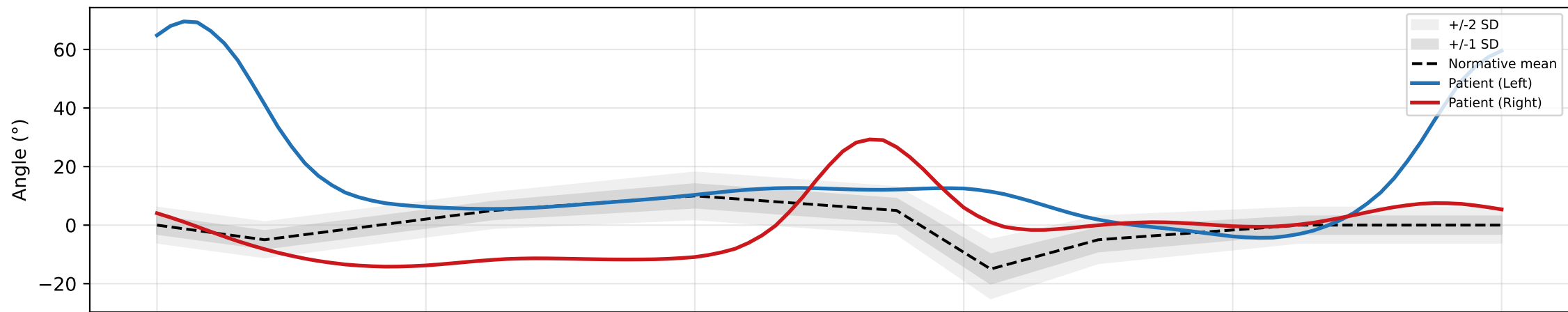

Trunk — vs Normative (adult)

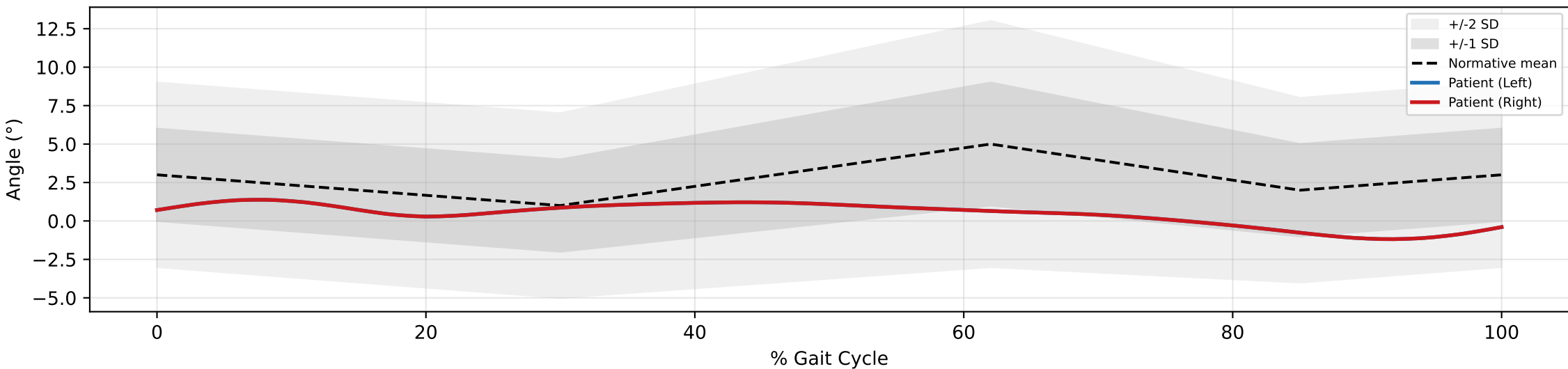


GVS / MAP Profile

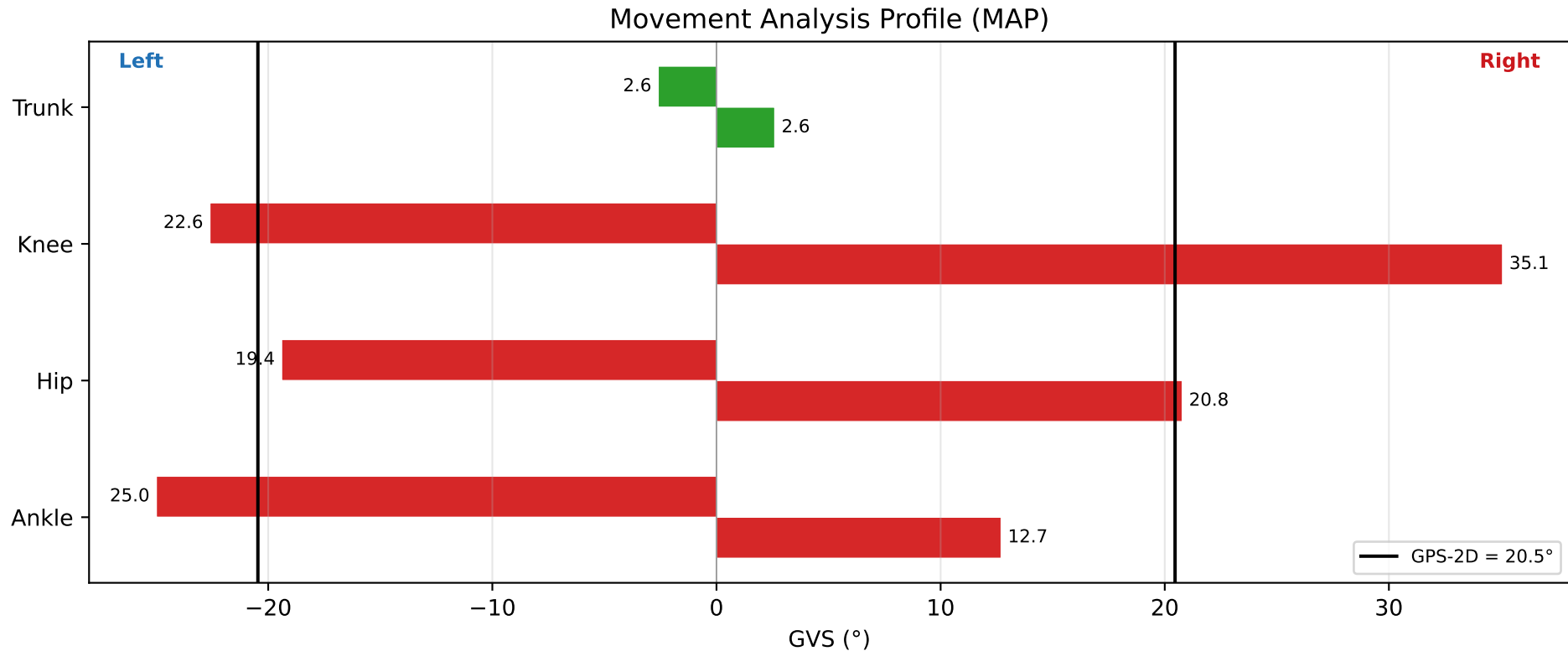


Quality Dashboard


Detection Rate per Landmark

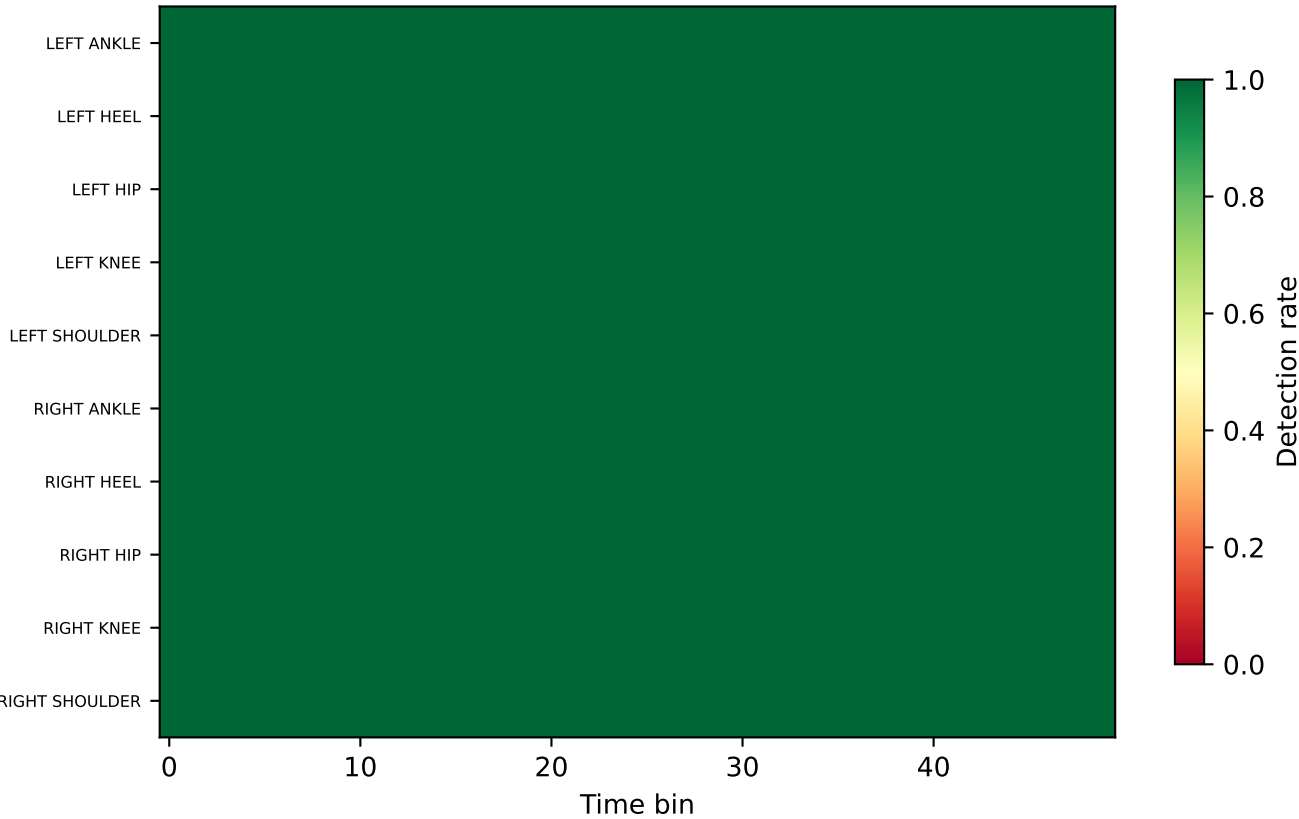

Per-frame Confidence

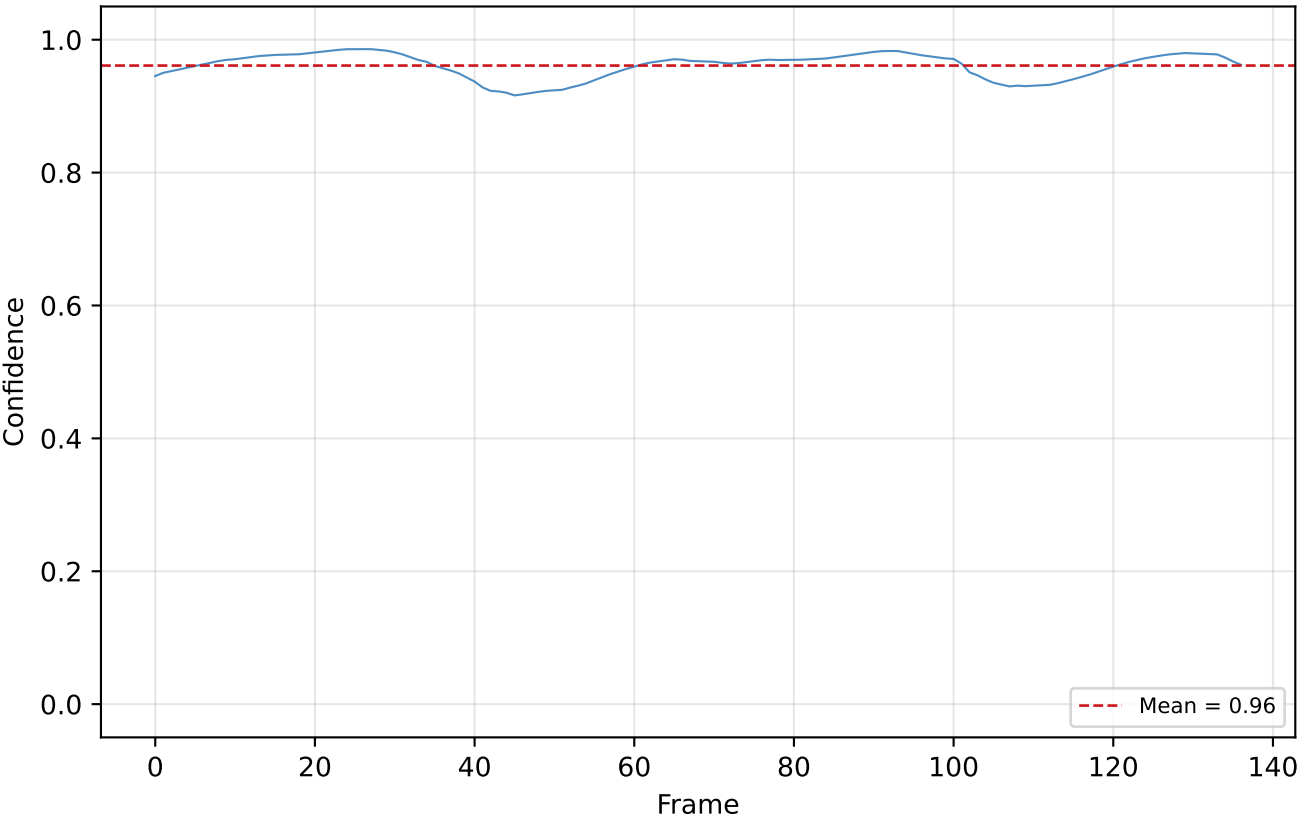

Gap Locations (0 gap frames)

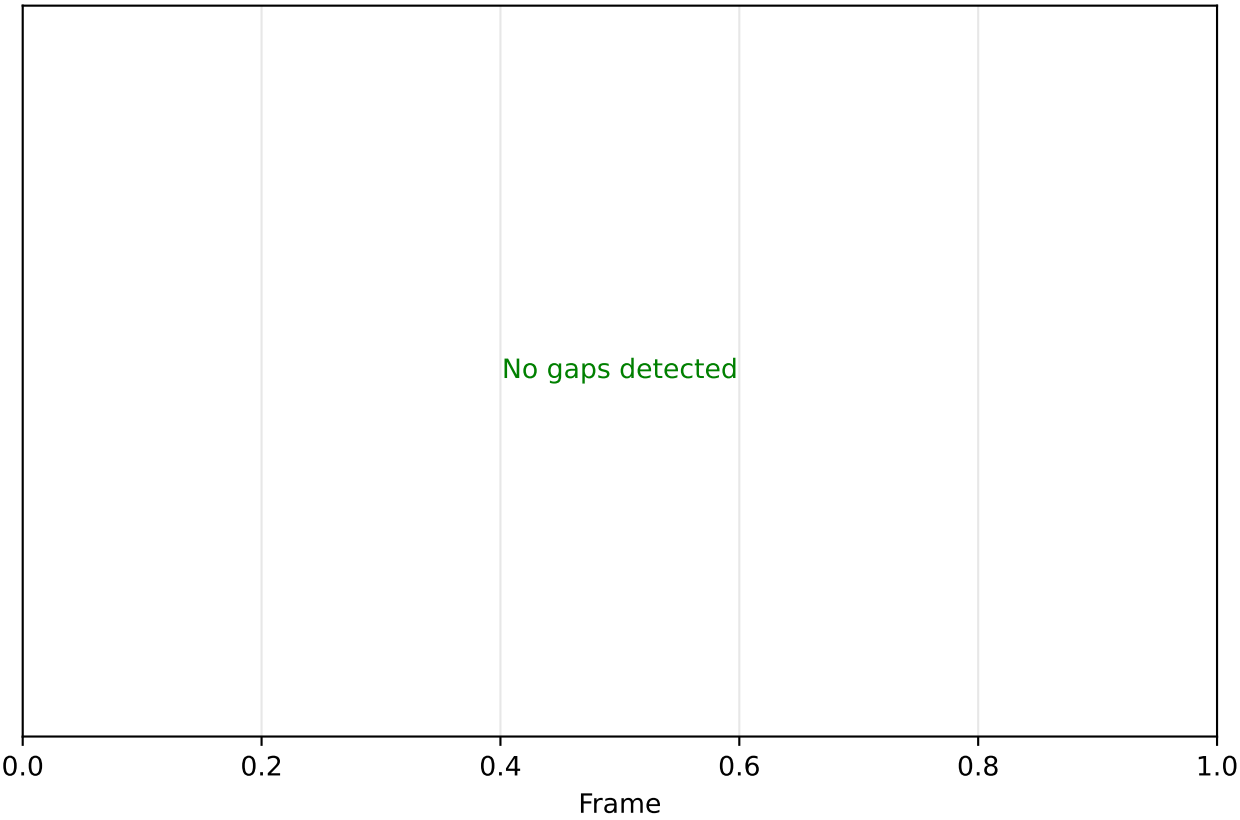

Overall Quality Score

99

Quality: GOOD

Detection rate: 100.0%
Mean confidence: 0.96
Gap percentage: 0.0%
Jitter score: 0.0008

Detailed Report

FULL REPORT – GAIT ANALYSIS

Source: /Users/adanmoallemi/Documents/MyoGait/videos/15-JUL-2026_pilot_ST_coronal_forward_1.MP4
FPS: 59.94005994005994 | Resolution: 1920x1080
Model: mediapipe | Correction: 0.8

SPATIO-TEMPORAL PARAMETERS

Cadence: 110.7 steps/min (normal: 100-120)
Stride (mean): 1.084 +/- 0.0 s
Step (mean): 0.542 +/- 0.292 s
Stance L: 53.8% Stance R: None% (normal: ~60%)
Double support: None% (normal: ~20%)

JOINT RANGES OF MOTION (ROM)

LEFT:
Hip 0.1 a 10.7 (ROM: 10.6)
Knee -9.1 a 31.2 (ROM: 40.3)
Ankle -4.3 a 69.6 (ROM: 73.9)
RIGHT:
Hip -3.7 a 7.6 (ROM: 11.3)
Knee -30.9 a -2.7 (ROM: 28.2)
Ankle -14.2 a 29.3 (ROM: 43.4)

SYMMETRY

hip_rom_si 6.2%
knee_rom_si 35.5%
ankle_rom_si 51.9%
step_time_si 0.0%
overall_si 33.6%
(SI < 10% = excellent, 10-20% = good, >20% = asymmetry)

VARIABILITY

Cycle duration CV: 0.0%
Stance % CV: 52.7%

ALERTS

* Asymmetry knee: SI=35.5% (>20%)
* Asymmetry ankle: SI=51.9% (>20%)
* Asymmetry stance_time: SI=74.5% (>20%)

Normal reference values:

Hip ROM: 40-50° | Knee ROM: 60-70° | Ankle ROM: 30-40°

Generated by myogait v0.6.0